
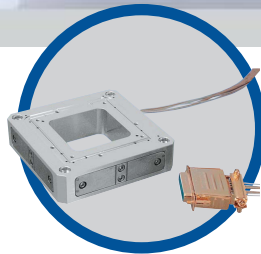


# Tools for Microscopy & Imaging

PDF Catalog, Updated July 2008

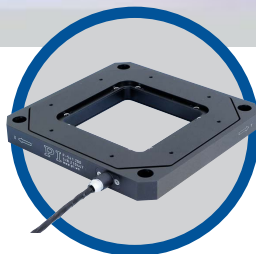

# Tools for Microscopy & Imaging

## Piezo Stages & Piezo Objective Scanners

### P-721.CLQ PIFOC® High-Speed Nanofocusing/Scanning Z-Drives



- Scans and Positions Objectives with Sub-nm Resolution
- High Linearity and Stability with Direct-Measuring Sensors
- Travel to 100  $\mu\text{m}$ , Fast Settling Time
- Very Straight Motion for Enhanced Focus Stability
- Ask about DIC Prism Holder Option
- Compatible with Metamorph™ Imaging Software
- Quick Lock Adapter for Easy Attachment

See p. 5 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

### P-725 PIFOC® Long-Range, High-Speed Nanofocusing Z-Drives



- Scans and Positions Objectives with Sub-nm Resolution
- High Linearity and Stability with Direct-Measuring Sensors
- Travel to 460  $\mu\text{m}$ , Fast Settling Time
- Frictionless Precision Flexure Guiding System
- Enhanced Guiding Precision for Better Focus Stability
- Ask about DIC Prism Holder Option
- Controller Compatible with Metamorph™ Imaging Software
- Quick Lock Adapter for Easy Attachment

See p. 9 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

### P-726 High-Power PIFOC® Nanofocusing Z-Drives



- Resonant Frequency: 430 Hz (with 300g), 1100 Hz Without Objective
- Closed-Loop Position Noise: <0.4 nm
- Position Stability Over 100s: 4nm
- Maximum Operating Frequency (amplitude <10  $\mu\text{m}$ ): 150 Hz
- Settling Time (10  $\mu\text{m}$  steps): 6 ms

See p. 11 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

### P-737 PIFOC® Z-Axis Microscopy Piezo Stage for High-Resolution Sample Positioning and Scanning




- High-Speed Piezo Motion with Travel Ranges to 500  $\mu\text{m}$
- Nanometer Resolution
- Large Clear Aperture to Accommodate Specimen Holders
- Perfect Mechanical Fit to OEM Manual or Motorized Stages
- Sub-Millisecond Response Times

See p. 12 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

# Tools for Microscopy & Imaging

## Piezo Stages & Piezo Objective Scanners


### M-686 PLine® XY Piezo Linear-Motor Stage

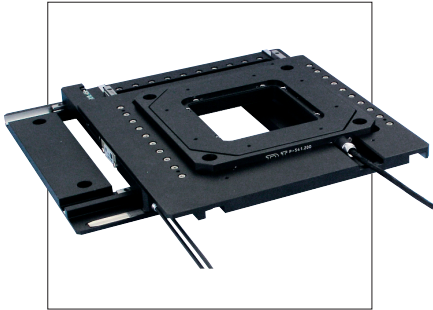

- Integrated Closed-Loop Piezomotor Drives Provide High Speed to 100 mm/s
- Travel Ranges 25 x 25 mm
- Integrated Linear Encoders with 0.1  $\mu\text{m}$  Resolution
- Compact Design: 32 mm Profile Height, 170 x 170 mm Footprint
- Clear Aperture 78 x 78 mm, 66 x 66 mm in Extreme Position
- Self-Locking at Rest
- Compatible with PI Piezo Nanopositioning / Scanning Stages

See p. 14 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

### P-541.2CD and P-733.2DD Low-Profile, XY Piezo Scanning Microscopy Stages

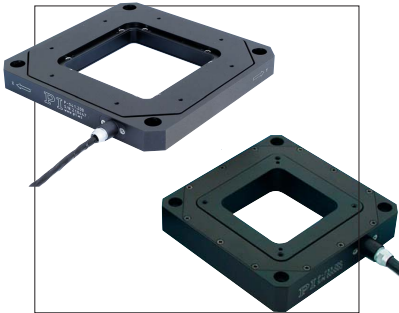

- Low Profile for Easy Integration
- Higher Speed Through Direct Drive
- Travel to 200 x 200  $\mu\text{m}$  or 30 x 30 High-Speed Version
- Frictionless Precision Flexure Guiding System
- Parallel-Kinematics/Metrology for Enhanced Responsiveness/Multi-Axis Precision
- 80 x 80 mm or 50 x 50 mm Clear Aperture
- PICMA® High-Performance Piezo Actuators for Superior Lifetime

See pp. 16-25, or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

### P-714 Low-Profile OEM XY Piezo-Scanners for Imaging Applications



- Ideal for Interlacing, CCD Resolution Enhancement
- Compact Size of only 45 x 45 x 6 mm
- Clear Aperture
- Highly Cost-Efficient Design
- 15 x 15  $\mu\text{m}$  Travel Range
- Parallel-Kinematics Design for Higher Dynamics

See p. 26 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

### P-733.2UD Non-Magnetic XY Scanning Stage for UHV to 10<sup>-9</sup> hPa



- Travel Ranges of 100 x 100  $\mu\text{m}$
- Resolution to 0.1 nm
- Direct Metrology with Capacitive Sensors
- Parallel Kinematics for Better Multi-Axis Accuracy and Dynamics
- Parallel-Motion Metrology for Highest Linearity and Stability
- 50 x 50 mm Clear Aperture for Transmitted-Light Applications

See p. 22 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

# Tools for Microscopy & Imaging

## Piezo Stages & Piezo Objective Scanners

### P-363 PicoCube® High-Speed, XY(Z) Piezo Stages for Nanotechnology, SPM, AFM



- Ultra-High-Performance Closed-Loop Scanner for AFM/SPM
- Compact Manipulation Tool for Nanotechnology
- Resonant Frequency 9.8 kHz
- Ultra-High-Precision Capacitive Feedback
- 50 Picometers Resolution, 5 x 5 x 5  $\mu\text{m}$  Travel Range
- Very Small Package, Rugged Design

See p. 28 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)


### PT120 - PT140 Piezoceramic Tubes (HVPZT) and P-151 PICA™-Shear Piezo Actuators



- Standard Custom Sizes
- For OEM Applications
- XYZ-Positioning
- Sub-Nanometer Resolution
- For Scanning Microscopy (STM, AFM, etc.)

See Online Datasheet Available at [www.pi.ws](http://www.pi.ws)


### E-625 and E-665 Piezo Amplifiers & Servo-Controllers

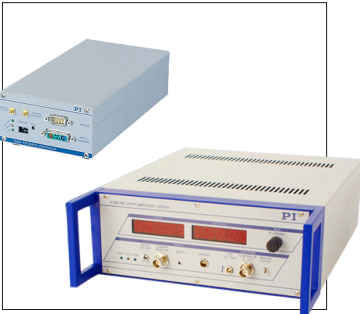

- Analog & 20-Bit High-Speed Digital Interface
- SGS, LVDT and Capacitive Sensor Servo-Control
- Wave Table for User-Defined Curves
- Network Capability with up to 4 or up to 12 Channels
- GCS (General Command Set) Compatible

See Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

### E-761 Digital Piezo Controller PCI Board for Piezo Stages

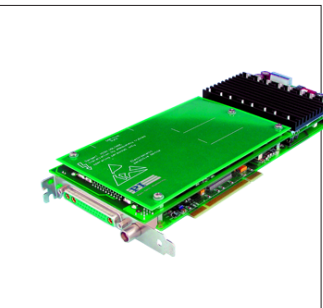

- For Piezo Stages with Capacitive Sensors
- High-Speed PCI Interface
- 3 Logical Axes, 4 Piezo Amplifiers
- Additional High-Bandwidth Analog Interface
- 32-Bit Digital Filters
- 24-Bit Ultra-Low-Noise DAC Converters
- Coordinate Transformation for Parallel-Kinematics / Parallel-Metrology Systems

See p. 30 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

# P-721

## PIFOC® High-Speed Nanofocusing/Scanning Z-Drives with Direct Metrology



P-721.CLQ with Quick Lock option P-721.12Q and objective (not included).

- Scans and Positions Objectives with Sub-nm Resolution
- High Linearity and Stability with Direct-Measuring Sensors
- Travel to 100 µm, fast Settling Time
- Very Straight Motion for Enhanced Focus Stability
- Ask about DIC Prism Holder Option
- Compatible with Metamorph™ Imaging Software
- Quick Lock Adapter for Easy Attachment

P-721 PIFOCs® are high-speed, piezo-driven microscope objective nanofocusing/scanning devices, providing a positioning and scanning range of 100 µm with sub-nanometer resolution and very high linearity. They were designed for tasks such as surface metrology or deconvolution microscopy (Z-stack acquisition). The frictionless, flexure guiding system provides enhanced precision for superior focus stability with fast response for rapid settling and scanning.

### Application Examples

- Scanning interferometry
- Surface structure analysis
- Disk drive testing
- Autofocus systems
- Confocal microscopy
- Biotechnology
- Semiconductor test equipment

### Superior Accuracy Through Direct-Motion Metrology with Capacitive Sensors

P-721.CLQ and .CDQ are equipped with direct-measuring capacitive position sensors. Unlike conventional sensors, capacitive sensors measure the actual distance between the fixed frame and the moving part of the stage. They detect errors contributed by all components in the drive train—from the actuator through the flexures to the platform. This results in higher motion linearity, long-term stability, phase fidelity, and—because external disturbances are seen by the sensor immediately—a stiffer, faster-responding servo-loop.

P-721.LLQ is a lower cost version, equipped with a direct-measuring LVDT sensor providing 10 nm resolution. A variety of analog and digital controllers is available to drive these units. The P-721.CDQ can be operated with the E-665

servo-controller (see page 6-30) through an analog or digital interface. The E-750 high-speed, digital NanoAutomation® controller, (p. 6-12) can also be used.

### Easy Installation

Attach the objective to the PIFOC®. The new Quick Lock thread adapters allow easy attachment to any of a variety of objective thread types.

Then screw a thread adapter onto the microscope and attach the PIFOC® with the Quick Lock system. Mounting does not require rotation of the PIFOC® unit, and the optical path is extended by only 12.5 mm (infinity corrected microscope required; extension tubes are available to adjust path lengths of other objectives on a turret). Thread adapters are ordered separately.

### Working Principle / Reliability

PIFOCs® are equipped with the award winning PICMA® piezo-drives, integrated into a

sophisticated flexure guiding system. The wire-EDM-cut flexures are FEA modeled for zero stiction, zero friction and exceptional guiding precision. The ceramic-encapsulated PICMA® drives are more robust than conventional piezo actuators, featuring superior lifetime and performance in both dynamic and static applications.


Because guidance, actuators and sensors are all frictionless and maintenance-free, these nanopositioning systems achieve outstanding levels of reliability.

### Other PIFOCs®

PIFOCs® are also available with up to 460 µm travel (P-725, see page 2-22) and for open-loop operation only (P-720, see page 2-20). Custom units for moving a whole turret are available on request.

### Notes

See the “Selection Guide” on p. 2-14 ff. for comparison with other nanopositioning systems.

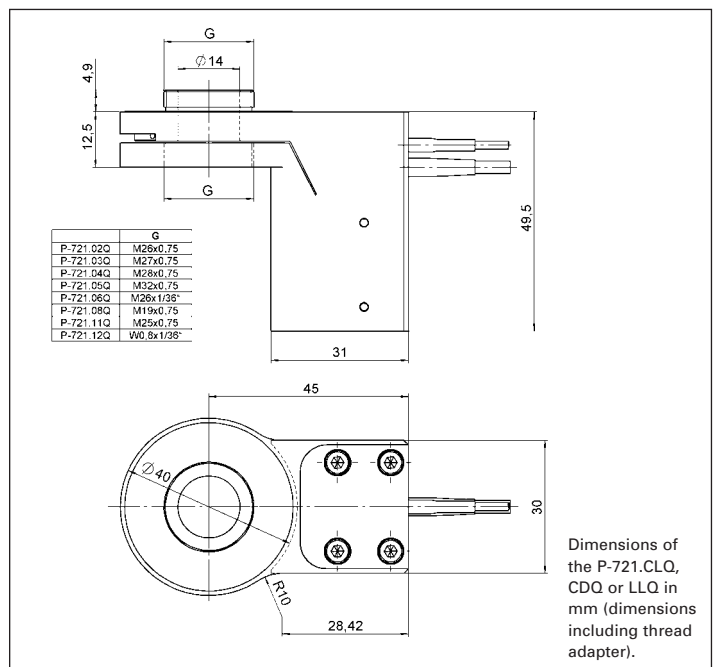

Dimensions of the P-721.CLQ, CDQ or LLQ in mm (dimensions including thread adapter).

## Ordering Information

### P-721.CLQ

PIFOC® Objective Positioner & Scanner, 100 µm, Capacitive Sensor, LEMO Connectors, for Quick Lock Thread Adapters

### P-721.CDQ

PIFOC® Objective Positioner & Scanner, 100 µm, Capacitive Sensor, Sub-D connector, for Quick Lock Thread Adapters

### P-721.LLQ

PIFOC® Objective Positioner & Scanner, 100 µm, LVDT Sensor, LEMO Connectors, for Quick Lock Thread Adapters

### P-721.0LQ

PIFOC® Objective Positioner & Scanner, 100 µm, Open-Loop, LEMO Connectors, for Quick Lock Thread Adapters

### Accessories: Quick Lock Adapters

#### P-721.11Q

Quick Lock Adapter, M25 x 0.75

#### P-721.12Q

Quick Lock Adapter, W0.8 x 1/36"

#### P-721.02Q

Quick Lock Adapter, M26 x 0.75

#### P-721.03Q

Quick Lock Adapter, M27 x 0.75

#### P-721.04Q

Quick Lock Adapter, M28 x 0.75

#### P-721.05Q

Quick Lock Adapter, M32 x 0.75

#### P-721.06Q

Quick Lock Adapter, M26 x 1/36"

#### P-721.08Q

Quick Lock Adapter, M19 x 0.75

### Extension Tubes for Objectives

#### P-721.90Q

12.5 mm Extens. Tube, W0.8 x 1/36"

#### P-721.91Q

12.5 mm Extens. Tube, M25 x 0.75

#### P-721.92Q

12.5 mm Extens. Tube, M26 x 0.75

#### P-721.93Q

12.5 mm Extens. Tube, M27 x 0.75

#### P-721.94Q

12.5 mm Extens. Tube, M28 x 0.75

#### P-721.95Q

12.5 mm Extens. Tube, M32 x 0.75


#### P-721.96Q

12.5 mm Extens. Tube, M26 x 1/36"

#### P-721.98Q

12.5 mm Extens. Tube, M19 x 0.75

Ask about custom designs!

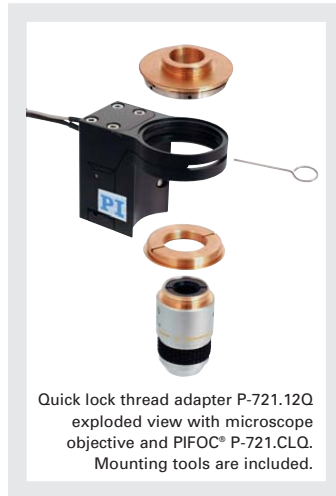

Quick lock thread adapter P-721.12Q exploded view with microscope objective and PIFO C P-721.CLQ. Mounting tools are included.



Custom PIFO C nosepiece nanopositioner for microscope turret.

## Technical Data

Models	P-721.CLQ P-721.CDQ	P-721.LLQ	P-721.0LQ	Units	Notes see p. 2-84
Max. objective diameter	39	39	39	mm	
Open-loop travel @ 0 to 100 V	90	90	90	µm ±20%	A2
Closed-loop travel	100	100	-	µm	A5
Integrated feedback sensor	Capacitive	LVDT	-		B
* Closed- / open loop resolution	0.7 / 0.5	10 / 0.5	- / 0.5	nm	C1
Closed-loop linearity (typ.)	0.03	0.1	-	%	
Full-range repeatability (typ.)	±5	±20	-	nm	C3
Stiffness	0.3	0.3	0.3	N/µm ±20%	D1
Push/pull force capacity (in operating direction)	100 / 20	100 / 20	100 / 20	N	D3
Tilt (θ <sub>X</sub> ) (typ.)	0.5	0.5	0.5	µrad	E1
Tilt (θ <sub>Y</sub> ) (typ.)	13	13	13	µrad	E1
Lateral runout (Y) (typ.)	100	100	100	nm	E2
Electrical capacitance	3.0	3.0	3.0	µF ±20%	F1
** Dynamic operating current coefficient (DOCC)	4.2	4.2	4.2	µA/(Hz x µm)	F2
Unloaded resonant frequency	580	580	580	Hz ±20%	G2
Resonant frequency @ 120 g load	250	250	250	Hz ±20%	G3
Resonant frequency @ 200 g load	190	190	190	Hz ±20%	G3
Operating temperature range	-20 to 80	-20 to 80	-20 to 80	°C	H2
Voltage connection	VL	VL	VL		J1
Sensor connection	2 x C	L	-		J2
Weight (with cables)	240	230	220	g ±5%	
Body material	Al	Al	Al		L
Recommended driver/ controller (codes explained p. 2-17)	H, M, F	H, E	G, C		

\* For calibration information see p. 2-8. Resolution of PI piezo nanopositioners is not limited by friction or stiction. The value given is noise equivalent motion with E-503 amplifier.

\*\* Dynamic Operating Current Coefficient in µA per Hz and µm. Example: Sinusoidal scan of 30 µm at 10 Hz requires approximately 1.2 mA drive current.

# P-720

## PIFOC® High-Speed Microscope Objective Nanofocusing/Scanning Z-Drives



P-720 objective nanofocusing / scanning drive (objective not included)

- Scans and Positions Objectives with Sub-nm Resolution
- Low Inertia for Fast Settling
- Frictionless Precision Flexure Guiding System
- Travel to 100  $\mu\text{m}$
- Straightness of Travel  $\leq 13 \mu\text{m}$
- PICMA® High-Performance Piezo Drives

P-720 PIFOCs® are high-speed, piezo-driven microscope objective nanofocusing/scanning devices which can be mounted on most microscopes. The frictionless, flexure guiding system combines high guiding precision for superior focus stability with fast response for rapid settling and scanning. The units are screwed between

the turret and the objective, providing a positioning and scanning range of up to 100  $\mu\text{m}$  with sub-nanometer resolution, while extending the optical path by only 13 mm (infinity-corrected microscope required; extension tubes are available to adjust path lengths of other objectives on the turret). The standard thread is W0.8 x 1/36", for alternate threads see the P-721.CLQ. For larger positioning ranges, to 460  $\mu\text{m}$ , see the P-725, page 2-22.

### Operation

The P-720 open-loop PIFOC® is designed for fast, high-resolution positioning and scanning tasks where the absolute position is not important or where an external sensor is used. The vertical position of the objective is roughly proportional to

the drive voltage (see p. 4-15 ff. in the "Tutorial: Piezoelectrics in Positioning" section for behavior of open-loop piezos). If absolute position control, high linearity, or repeatability in the nanometer range is required, refer to the P-721 and P-725 closed-loop devices (see pages 2-20 and 2-22).

### Working Principle / Reliability

PIFOCs® are equipped with the award winning PICMA® piezo drives, integrated into a sophisticated flexure guiding system. The wire-EDM-cut flexures are FEA modeled for zero stiction, zero friction and exceptional guiding precision. The ceramic-encapsulated PICMA® drives are more robust than conventional piezo actuators, featuring superior lifetime and performance in both dynamic and static applications.

### Notes

See the "Selection Guide" on p. 2-14 ff. for comparison with other nanopositioning systems.

### Ordering Information

**P-720.00**  
PIFOC® Objective Positioner & Scanner, 100  $\mu\text{m}$ , W0.8 x 1/36"


For PIFOC® Objective Positioners & Scanners with direct metrology and travel ranges to 400  $\mu\text{m}$  see P-721 and P-725, p. 2-20 and p. 2-22

**P-720.01**  
Objective extension tube, 13 mm


Ask about custom designs!

### Application Examples

- Scanning interferometry
- Surface structure analysis
- Disk drive testing
- Autofocus systems
- Confocal microscopy
- Biotechnology
- Semiconductor test equipment

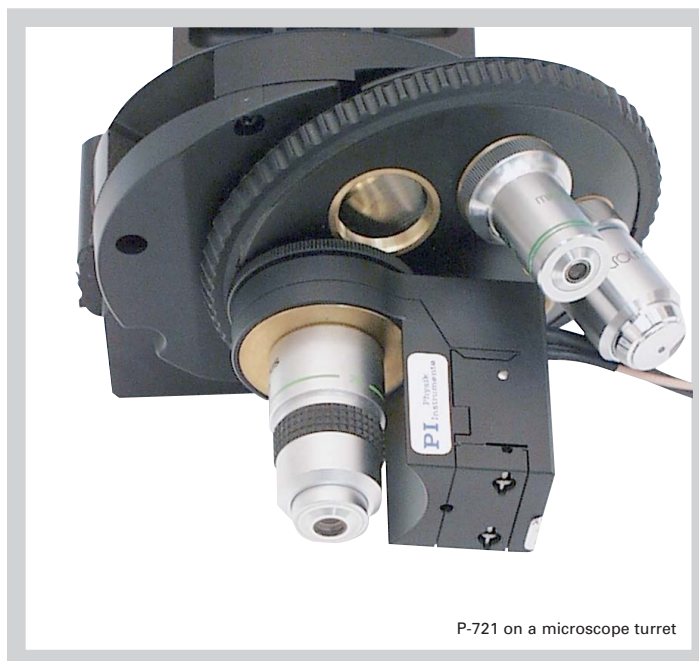

P-721 on a microscope turret

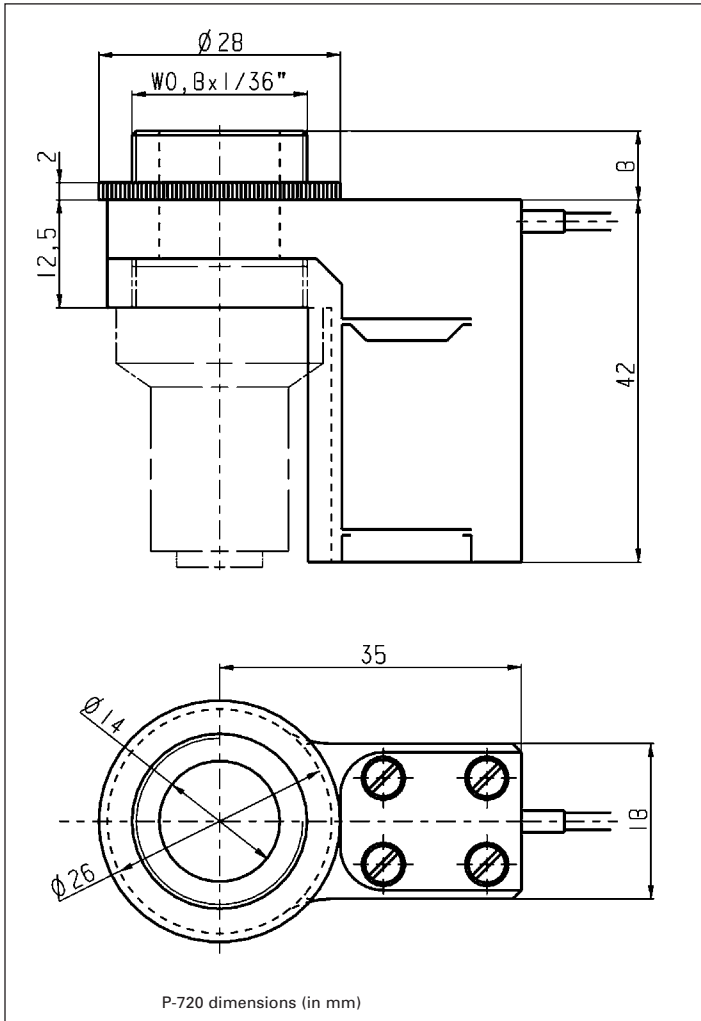

### Technical Data

Models	P-720.00	Units	Notes see p. 2-84
Max. objective diameter	25	mm	
Open-loop travel @ 0 to 100 V	100	$\mu\text{m} \pm 20\%$	A2
* Open-loop resolution	1	nm	C1
Stiffness	0.2	$\text{N}/\mu\text{m} \pm 20\%$	D1
Push/pull force capacity (in operating direction)	100 / 20	N	D3
Tilt ( $\theta_x, \theta_y$ ) (typ.)	13	$\mu\text{rad}$	E1
Lateral runout (Y) (typ.)	100	nm	E2
Electrical capacitance	3.0	$\mu\text{F} \pm 20\%$	F1
** Dynamic operating current coefficient (DOCC)	3.8	$\mu\text{A}/(\text{Hz } 3 \mu\text{m})$	F2
Unloaded resonant frequency	400	$\text{Hz} \pm 20\%$	G2
Resonant frequency @ 120 g load	180	$\text{Hz} \pm 20\%$	G3
Resonant frequency @ 200 g load	150	$\text{Hz} \pm 20\%$	G3
Operating temperature range	-20 to 80	$^{\circ}\text{C}$	H2
Voltage connection	VL		J1
Weight (with cables)	100	$\text{g} \pm 5\%$	
Body material	Al		L
Recommended driver/controller (codes explained p. 2-17)	G, C, (A)		


\* For calibration information see p. 2-8.  
Resolution of PI piezo nanostages is not limited by friction or stiction. The value given is noise equivalent motion with E-503 amplifier.

\*\* Dynamic Operating Current Coefficient in  $\mu\text{A}$  per Hz and  $\mu\text{m}$ . Example: Sinusoidal scan of  $30 \mu\text{m}$  at 10 Hz requires approximately 1.2 mA drive current.



# P-725

## PIFOC® Long-Range, High-Speed Nanofocusing Z-Drives with Direct Metrology

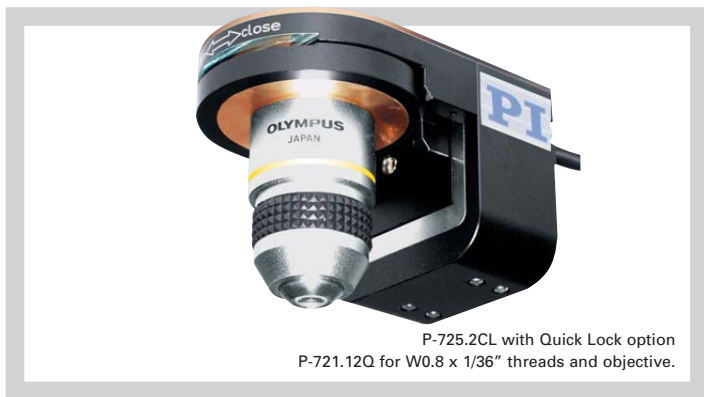

P-725.2CL with Quick Lock option  
P-721.12Q for W0.8 x 1/36" threads and objective.

- Scans and Positions Objectives with Sub-nm Resolution
- High Linearity and Stability with Direct-Measuring Capacitive Sensors
- Travel to 460 μm, Fast Response & Settling Time
- Frictionless Precision Flexure Guiding System
- Enhanced Guiding Precision for Better Focus Stability
- Ask about DIC Prism Holder Option
- Controller Compatible with Metamorph™ Imaging Software
- Quick Lock Adapter for Easy Attachment

P-725 PIFOCs® are long-travel, high-speed, piezo-driven microscope objective nanofocusing/scanning devices. Despite the increased travel ranges (up to 460 μm), they are 20% shorter than P-721 units (see page 2-20), while providing sub-nanometer resolution. Equipped with ultra-precise, direct-measuring capacitive sensors, these devices are ideal for tasks such as surface metrology or deconvolution microscopy (Z-stack acquisition). The newly designed, frictionless, flexure guiding system provides enhanced preci-

sion for superior focus stability with fast response for rapid settling and scanning.

P-725 PIFOCs® are mounted between the turret and the objective, extending the optical path by only 12.5 mm (infinitely corrected microscope required; extension tubes are available to adjust path lengths of other objectives on the turret). Custom designs for positioning the complete turret are available on request.

### Superior Accuracy Through Direct-Motion Metrology with Capacitive Sensors

P-725s are equipped with absolute-measuring, direct-metrology capacitive sensors. These sensors make possible motion linearity to 0.03% with effective resolution in the sub-nanometer range. They boast high bandwidth and exhibit no periodic errors. Unlike conventional sensors, capacitive sensors measure the actual distance between the fixed frame and the mov-

ing part of the stage. They detect errors contributed by all components in the drive train—from the actuator through the flexures to the platform. This results in higher motion linearity, long-term stability, phase fidelity, and—because external disturbances are seen by the sensor immediately—a stiffer, faster-responding servo-loop. See p. 2-4 *ff.* and p. 5-2 *ff.* for more information.

### Fastest Step-and-Settle: 25 Milliseconds for 250 Microns

The P-725.2CL can perform a 250 μm step to 1% accuracy in only 25 ms (E-665.CR controller, no load) and in 50 ms with a load of 150 g.

### Controllers

A variety of analog and digital controllers is available to drive these units. The P-725.xCD models can be operated with the E-665 servo-controller through an analog or digital interface (see page 6-30). The E-750 high-speed digital NanoAutomation® controller can also be used (see page 6-12).

### Working Principle / Reliability

PIFOCs® are equipped with the award winning PICMA® piezo drives, integrated into a sophisticated flexure guiding system. The wire-EDM-cut flexures are FEA modeled for zero stiction, zero friction and exceptional guiding precision. The ceramic-encapsulated

### Ordering Information

- P-725.1CL\***  
PIFOC® Objective Positioner & Scanner, 100 μm, Capacitive Sensor, LEMO Connectors, for Quick Lock Thread Adapters
  - P-725.1CD**  
PIFOC® Objective Positioner & Scanner, 100 μm, Capacitive Sensor, Sub-D Connector, for Quick Lock Thread Adapters
  - P-725.2CL\***  
PIFOC® Objective Positioner & Scanner, 250 μm, Capacitive Sensor, LEMO Connectors, for Quick Lock Thread Adapters
  - P-725.2CD**  
PIFOC® Objective Positioner & Scanner, 250 μm, Capacitive Sensor, Sub-D Connector, for Quick Lock Thread Adapters
  - P-725.4CL\***  
PIFOC® Objective Positioner & Scanner, 400 μm, Capacitive Sensor, LEMO Connectors, for Quick Lock Thread Adapters
  - P-725.4CD**  
PIFOC® Objective Positioner & Scanner, 400 μm, Capacitive Sensor, Sub-D Connector, for Quick Lock Thread Adapters
- \*Also available w/o sensor (open-loop): P-725.10L, P-725.20L and P-725.40L.

**Accessories**  
Quick Lock adapters and objective distance tubes see p. 2-21


**Ask about custom designs!**

PICMA® drives are more robust than conventional piezo actuators, featuring superior lifetime and performance in both dynamic and static applications.

Because guidance, actuators and sensors are all frictionless and maintenance-free, these

### Application Examples

- Scanning interferometry
- Surface structure analysis
- Disk drive testing
- Autofocus systems
- Confocal microscopy
- Biotechnology
- Semiconductor test equipment

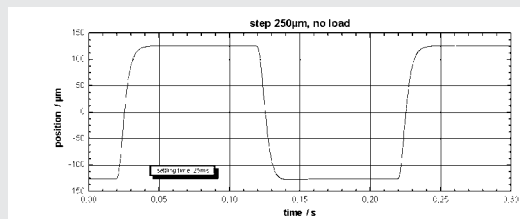



25 ms for a 250 μm step top dynamic performance of the P-725.2CL

nanopositioning systems achieve outstanding levels of reliability.

### Notes

For low-profile, large aperture Z stages, see the P-541.ZSL on p. 2-48 and the P-541.ZCD on p. 2-50.

See the “Selection Guide” on p. 2-14 ff. for comparison with other nanopositioning systems.

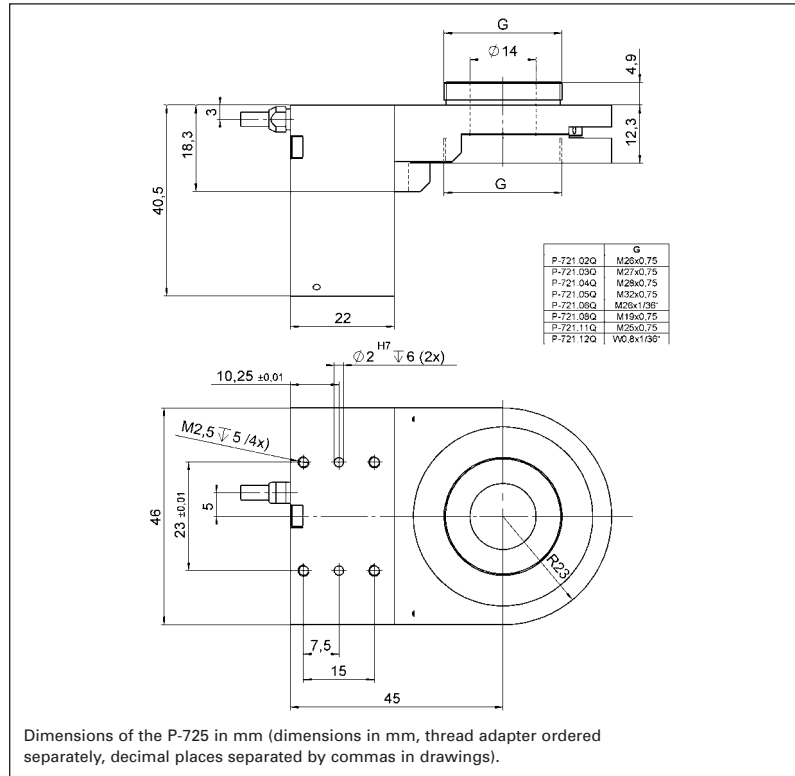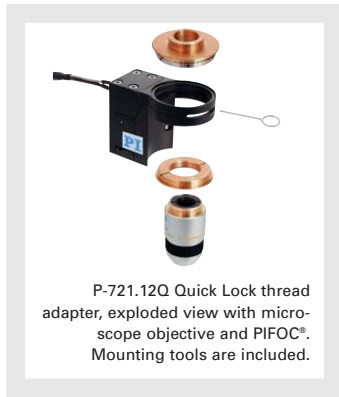

### Technical Data

Models	P-725.1CL, P-725.1CD	P-725.2CL, P-725.2CD	P-725.4CL, P-725.4CD	Units	Notes see p. 2-84
Max. objective diameter	39	39	39	mm	
Min. open-loop travel @ -20 to 120 V	150	330	460	$\mu\text{m} \pm 20\%$	A2
Closed-loop travel	100	250	400	$\mu\text{m}$	A5
Integrated feedback sensor	Capacitive	Capacitive	Capacitive		B
* Closed-loop resolution	0.65	0.75	1.25	nm	C1
Closed-loop linearity (typ.)	0.03	0.03	0.03	%	
Full-range repeatability (typ.)	$\pm 2$	$\pm 3$	$\pm 5$	nm	C3
Stiffness	0.25	0.2	0.15	$\text{N}/\mu\text{m} \pm 20\%$	D1
Push/pull force capacity (in operating direction)	100 / 20	100 / 20	100 / 20	N	D3
Tilt ( $\theta_x$ ) (typ.)	1	6	18	$\mu\text{rad}$	E1
Tilt ( $\theta_y$ ) (typ.)	20	45	25	$\mu\text{rad}$	E1
Lateral runout (Y) (typ.)	40	40	40	nm	E2
Electrical capacitance	4.2	6.0	6.0	$\mu\text{F} \pm 20\%$	F1
** Dynamic operating current coefficient (DOCC)	4.0	2.5	1.9	$\mu\text{A}/(\text{Hz} \times \mu\text{m})$	F2
Unloaded resonant frequency	530	330	200	$\text{Hz} \pm 20\%$	G2
Resonant frequency @ 120 g load	205	180	115	$\text{Hz} \pm 20\%$	G3
Resonant frequency @ 200 g load	160	140		$\text{Hz} \pm 20\%$	G3
Operating temperature range	-20 to 80	-20 to 80	-20 to 80	$^{\circ}\text{C}$	H2
Voltage connection	VL ***	VL ***	VL ***		J1
Sensor connection	2 x C ***	2 x C ***	2 x C ***		J2
Weight (with cables)	215 (230)	230 (245)	230 (245)	$\text{g} \pm 5\%$	
Body material	Al	Al	Al		L
Recommended driver/controller (codes explained p. 2-17)	H, M, F	H, M, F	H, M, F		

\* For calibration information see p. 2-8.  
Resolution of PI piezo nanopositioners is not limited by friction or stiction. The value given is noise equivalent motion with E-503 amplifier.

\*\* Dynamic Operating Current Coefficient in  $\mu\text{A}$  per Hz and  $\mu\text{m}$ . Example (P-725.2C): Sinusoidal scan of  $30 \mu\text{m}$  at 10 Hz requires approximately 0.75 mA drive current.

\*\*\* P-725.xCD with Sub-D connector for voltage and sensor.

## News & Applications

# The fastest microscope objective Z-piezo stage P-726 High-Performance PIFOC®




PI has extended its successful PIFOC® line of objective nanopositioners with a new, high-power unit. The P-726 PIFOC® can move heavy, high-NA objectives (>200 g, 7 oz) faster than any other system currently available, with accuracies of 1 nm over a travel range of 100 µm.

Here, the convincing performance specifications – achieved with a load of 300 g (11 oz): The E-726 can be operated with a variety of controllers ranging from the economical E-621 card, the E-625 controller to a number of digital controllers for the highest possible performance.

- Resonant frequency: 430 Hz (with 300 g), 1100 Hz without objective
- Closed-loop position noise: <0.4 nm
- Position stability over 100 s: 4 nm
- Maximum operating frequency (amplitude <10 µm): 150 Hz\*
- Maximum operating frequency with full travel (100 µm): 60 Hz\*
- Settling time (10 µm steps): 6 ms\*\*
- Settling time (full travel): 15 ms\*\*

\* Amplifier power limited.  
Higher performance with E-505,  
200 W power module  
\*\* Faster settling with smaller load

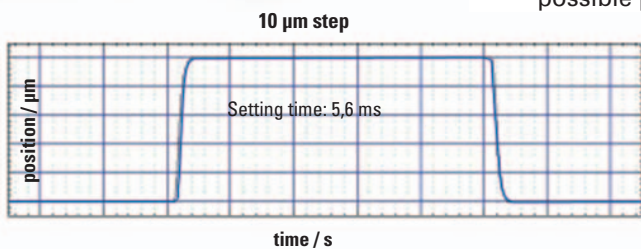

Settling behavior of the P-726 for a 10 µm step.

# P-737

## PIFOC® Z-Axis Microscopy Piezo Stage for High-Resolution Sample Positioning and Scanning



P-737 piezo Z-stage for high-resolution microscopy, shown with multiwell plate

- High-Speed Piezo Z-Motion with Travel Ranges to 250 µm
- Nanometer Resolution
- Large Clear Aperture to Accommodate Specimen Holders
- Perfect Mechanical Fit to XY OEM Manual or Motorized Stages
- Sub-Millisecond Response Times

PIFOC® P-737 high-speed vertical positioning systems are designed for use with XY microscopy stages – OEM manual stages as well as aftermarket motorized stages.

While the XY stage positions the sample, the piezo-actuator-based P-737 moves the sample along the optical axis to quickly and precisely adjust the focus. Vertical stepping with an accuracy in the nanometer range takes only a few milliseconds.

### Application Examples

- Fluorescence microscopy
- Confocal microscopy
- Biotechnology
- 3D Imaging
- Autofocus systems
- Medical technology

The large aperture is designed to accommodate a variety of specimen holders including slides or multiwell plates.

### High-Speed Z-Steps for Fast Focus Control and Z-Stack Acquisition

The immediate response of the solid-state piezo drives enables rapid Z-steps with typically 10 to 20 times faster step-and-settle times than classical stepper motor drives. This leads to more image acquisition speed and throughput.

### Closed-Loop Position Control for High-Precision and Stability

For high stability and repeatability, P-737 stages are equipped with position sensors. High-resolution, absolute measuring strain gauge sensors (SGS) are applied to appropriate places on the drive train and feed the platform

### Ordering Information

**P-737.1SL**  
PIFOC® Nanofocusing Z-Stage for Microscope Sample Holder, 100 µm, SGS, LEMO Connector, for Märzhäuser Microscope Stages

**P-737.2SL**  
PIFOC® Nanofocusing Z-Stage for Microscope Sample Holder, 200 µm, SGS, LEMO Connector, for Märzhäuser Microscope Stages


Versions with up to 500 µm travel or with direct-measuring, high-resolution capacitive sensors on request.

**Ask about custom designs!**

position information back to a piezoelectric controller. The sensors are connected in a full-bridge configuration to eliminate thermal drift, and assure optimal position stability and rapid response with nanometer resolution.

### Excellent Guiding Accuracy

Flexures optimized with Finite Element Analysis (FEA) are

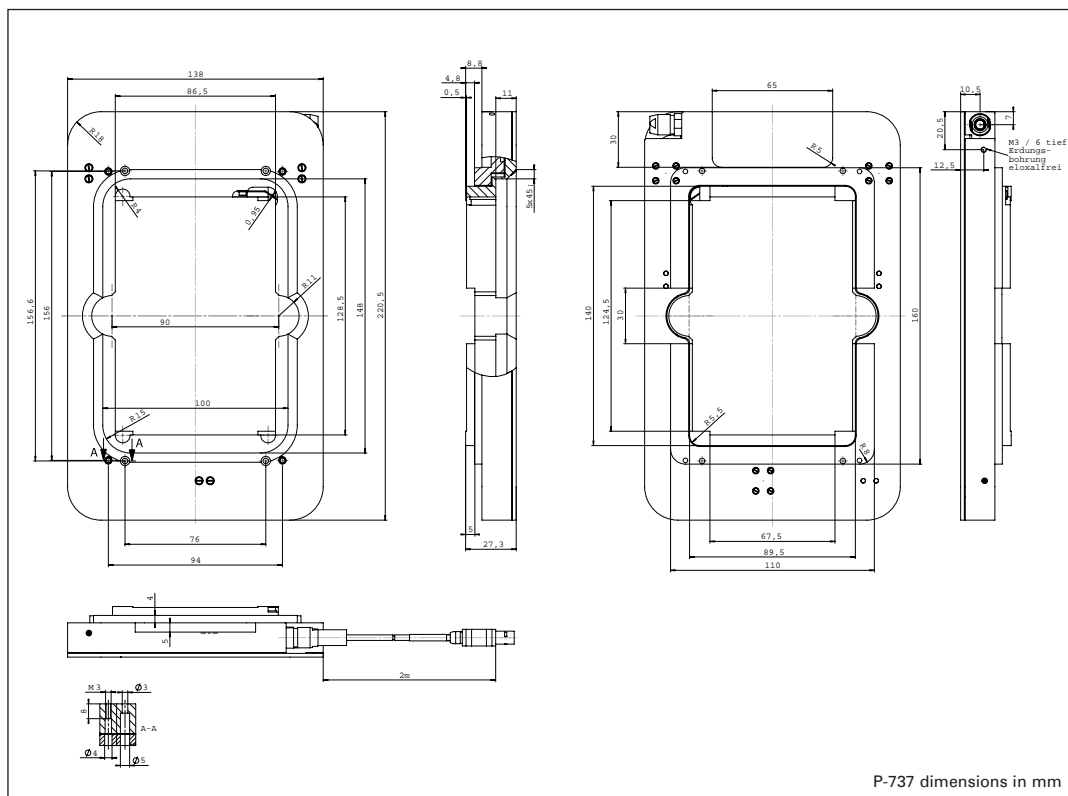


used to guide the stage. The FEA techniques give the design the highest possible stiffness in, and perpendicular to, the direction of motion, and minimize linear and angular runout. Flexures allow extremely high-precision motion, no matter how minute, as they are completely free of play and friction.

### Superior Lifetime with Ceramic-Encapsulated Piezos

Highest possible reliability is assured by the use of award-winning PICMA® multilayer piezo actuators. PICMA® actua-

tors are the only ceramic-encapsulated PZT actuators on the market, which makes them resistant to ambient humidity and leakage-current failures. They are thus far superior to conventional actuators in reliability and lifetime.

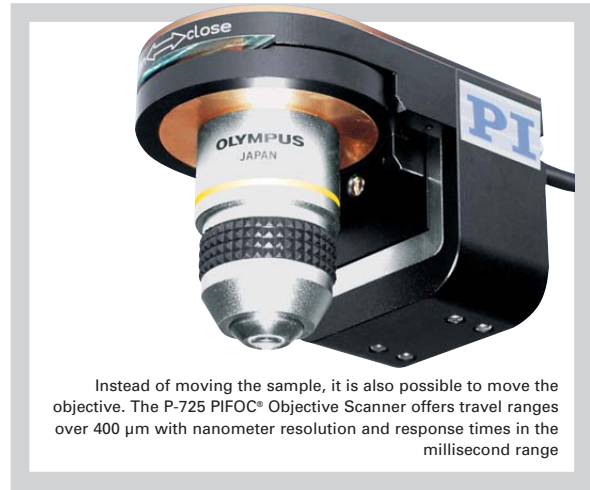

Instead of moving the sample, it is also possible to move the objective. The P-725 PIFOC® Objective Scanner offers travel ranges over 400  $\mu\text{m}$  with nanometer resolution and response times in the millisecond range

### Technical Data

Models	P-737.1SL	P-737.2SL	Units	Tolerance
Active axes	Z	Z		
<b>Motion and positioning</b>				
Integrated sensor	SGS	SGS		
Closed-loop travel	100	250	$\mu\text{m}$	calibrated
Closed-loop resolution	2.5	4	nm	typ.
Open-loop resolution	0.3	0.5	nm	typ.
Closed-loop linearity	0.2	0.2	%	typ.
Repeatability	6	12	nm	typ.
<b>Mechanical properties</b>				
Unloaded resonant frequency	400	250	Hz	$\pm 20\%$
Resonant frequency under load	350 (100 g)	210 (100 g)	Hz	$\pm 20\%$
Push/pull force capacity in motion direction	100 / 20	100 / 20	N	max.
Load	20	20	N	max.
<b>Drive properties</b>				
Piezo ceramic type	PICMA®	PICMA®		
Electrical capacitance	4.5	9.3	$\mu\text{F}$	$\pm 20\%$
Dynamic operating current coefficient	6.2	5.8	$\mu\text{A}/(\text{Hz} \times \mu\text{m})$	$\pm 20\%$
<b>Miscellaneous</b>				
Operating temperature range	-20 to 80	-20 to 80	$^{\circ}\text{C}$	
Material	Aluminum	Aluminum		
Dimensions	220.5 x 138 x 27.3	220.5 x 138 x 27.3	mm	
Mass	0.7	0.7	kg	$\pm 5\%$
Cable length	2.0	2.0	m	$\pm 10\text{ mm}$
Voltage connection	LEMO connector (low voltage)	LEMO connector (low voltage)		
Sensor connection	LEMO connector (SGS sensor)	LEMO connector (SGS sensor)		
Recommended controller/amplifier	E-625.SR, 14 W, single-channel desktop unit E-665.SR, 36 W, single-channel desktop unit E-610.SR, 14 W, OEM electronics E-500 modular system with output to 100 W			
<b>System properties</b>				
System configuration	E-500 system with E-503.00 14 W amplifier and E-509.S1 servo-control module for SGS			
Small-signal bandwidth	60	45	Hz	typ.
Settling time (to within 10% of step size)	10	15	ms	typ.

## M-686 PLine® XY Piezo Linear-Motor Stage

### Fast, Low Profile and Large Aperture with Direct Position Measurement



The M-686.D64 open-frame stage with closed-loop piezo motors provides 25 x 25 mm travel range

- **Integrated Closed-Loop Piezomotor Drives Provide High Speed to 100 mm/s**
- **Travel Ranges 25 x 25 mm**
- **Integrated Linear Encoders with 0.1 µm Resolution**
- **Compact Design:**  
32 mm Profile Height, 170 x 170 mm Footprint
- **Clear Aperture 78 x 78 mm, 66 x 66 mm in Extreme Position**
- **Self-Locking at Rest**
- **Compatible with PI Piezo Nanopositioning / Scanning Stages**

M-686 open-frame piezomotor stages are mainly designed for automated positioning applications in microscopy. The optimized form factor with a low profile height of only 32 mm and the standardized mounting pattern allows the combination with many PI standard nanopositioning systems.

#### Application Examples

- Biotechnology
- Microscopy
- Scanning microscopy
- Confocal microscopy
- Semiconductor testing
- Handling

#### Space Saving Piezomotors

Compared to conventional motorized translation stages, the M-686 provides a lower profile and smaller footprint. The compact PLine® piezoelectric linear motors and high-resolution linear encoders make both, the lead screw duct and the flanged, bulky stepper motor employed in traditional stages obsolete. In addition, the piezomotors are self-locking at rest and hold the stage in a stable position without heating up.

#### Compatibility to PI Nanopositioning and Scanning Stages

A number of standard PI piezo flexure stages (150 x 150 mm footprint) can be mounted directly on the M-686 open-

frame stage. Depending on the application, these highly specialized, ultra-precise nanopositioning systems are available as fast XY scanners (for fluorescence microscopy), as vertical Z positioners (3D imaging), or with up to 6 degrees of freedom.

#### Limit and Reference Switches

For the protection of your equipment, non-contact Hall-effect limit and reference switches are installed. The direction-sensing reference switch supports advanced automation applications with high precision.

#### Advantages of PLine® Micropositioning Systems

The ultrasonic piezoceramic drives used in PLine® micropositioners have a number of advantages over classical drives:

- Higher Accelerations, up to 5 g
- Speeds up to 500 mm/s
- Small Form Factor
- Self-Locking When Powered Down
- No Shafts, Gears or Other Rotating Parts
- Non-Magnetic and Vacuum-Compatible Drive Principle

#### Ordering Information

**M-686.D64**  
XY Open-Frame Stage with Closed-Loop PLine® Piezomotor Drives, 25 x 25 mm, 7 N, 0.1 µm Linear Encoder

**Ask about custom designs!**

#### Notes

Nanopositioning stages that fit directly on the M-686:

#### P-561 to P-563

PIMars™ XYZ Nanopositioning systems with up to 300 µm travel

#### P-541.2 to P-542.2


Low-profile microscopy XY scanners

#### P-541.Z

Low-profile Z/tip/tilt piezo nanopositioning stages for microscopy



Customized M-686 stage with a bigger footprint, to sink the piezo scanner by 10 mm. The system height together with the P-541 piezo scanner is reduced to only 33 mm

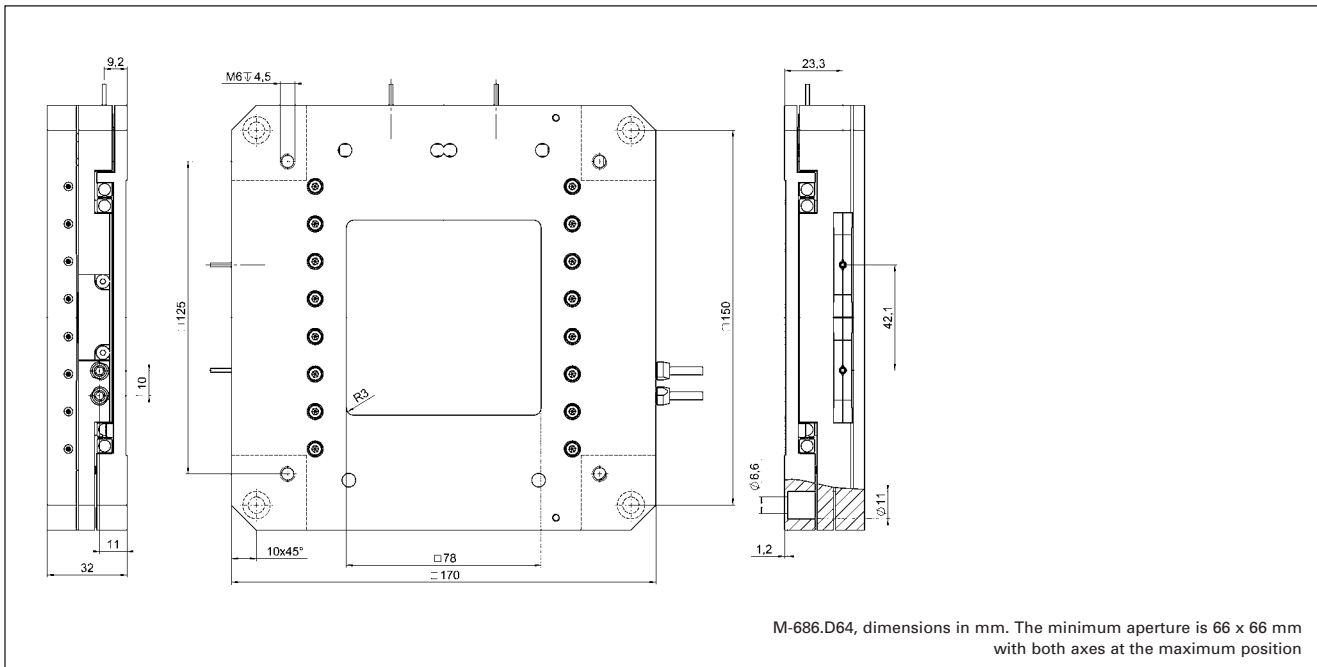

### Technical Data

	<b>M-686.D64</b>
Active axes	XY
<b>Motion and positioning</b>	
Travel range	25 x 25 mm
Integrated sensor	Linear encoder
Sensor resolution	0.1 $\mu\text{m}$
Design resolution	0.1 $\mu\text{m}$
Min. incremental motion	0.3 $\mu\text{m}$
Bidirectional repeatability	0.3 $\mu\text{m}$
Pitch / yaw	$\pm 50 \mu\text{rad}$
Max. velocity	100 mm/s
<b>Mechanical properties</b>	
Load Capacity*	50 N
Max. push/pull force	7 N
Max. lateral force	4 N
<b>Drive properties</b>	
Motor type	2 x PLine® P-664 per axis
Operating voltage	190 V (Peak-Peak)**
	67 V (RMS)**
Electrical power	10 W / axis***
<b>Miscellaneous</b>	
Operating temperature range	-20 to +50 °C
Material	Aluminium (black anodized)
Mass	1.2 kg
Cable length	1.5 m
Connector	2 x MDR connector, 14-pin
Recommended controller/driver	2 x C-866.D64 single-axis controller / driver 2 x C-185.D64 single-axis drive electronics


\* 10 N for max. velocity

\*\* The operating voltage or the piezomotor is supplied by the drive electronics which requires 12 VDC

\*\*\* For drive electronics



M-686 open-frame stage with P-541.2DD piezo scanner on top, providing a resolution of 0.1 nm and a scanning range of 30 x 30  $\mu\text{m}$ . The system height of the combination with the P-541 XY (or Z) piezo scanner is only 48 mm

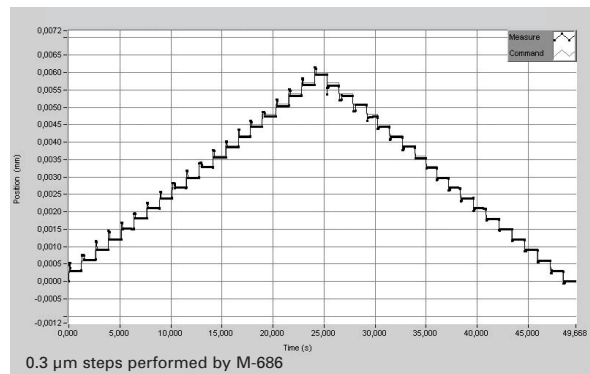

0.3  $\mu\text{m}$  steps performed by M-686

# P-541.ZCD

## Low-Profile Z/Tip/Tilt Piezo Nanopositioning Stages for Microscopy with Parallel Metrology



P-541 nanopositioning Z- and Z-tip/tilt stages feature a very low profile of 16.5 mm, a large 80 x 80 mm aperture and deliver highly accurate motion with sub-nanometer resolution. Strain gauge and capacitive sensors are available for feedback.

- Lower Profile for Easy Integration: 16.5 mm
- Z and Tip/Tilt Versions
- Choice of Strain Gauge (Lower Cost) or Capacitive Sensors (Higher Performance)
- 100 µm Linear Travel Range, 1 mrad Tilt
- 80 x 80 mm Clear Aperture
- PICMA® High-Performance Piezo Actuators for Superior Lifetime
- Ideal for fast Z-Stack Acquisition and 3D Deconvolution Imaging

### Low Profile, Optimized for Microscopy Applications

P-541 series nanopositioning Z-stages and Z-tip/tilt stages are designed for alignment, nano-focusing or metrology tasks, etc. They feature a very low profile of 16.5 mm, a large 80 x 80 mm aperture, and offer highly accurate motion with sub-nanometer resolution. XY

#### Application Examples

- Scanning microscopy
- Mask & wafer alignment
- Scanning interferometry
- Surface metrology
- Biotechnology
- Micromanipulation

scanning stages with the same footprint are also available (see p. 2-62).

Due to the single-module, low-profile design, the stages can easily be integrated in high-resolution microscopes.

### Superior Accuracy Through Direct-Motion Metrology with Capacitive Feedback Sensors

Integrated capacitive sensors measure the actual distance between the fixed frame and the moving part of the stage directly and include any flex or other errors in the drive train—from the actuator through the lever and flexures to the platform—in the measurement. This results in higher motion linearity, long-term stability,

phase fidelity, and—because external disturbances are seen by the sensor immediately—a stiffer, faster-responding servo-loop. See p. 2-4 *ff.* and p. 5-2 for more information.

### Dynamic Digital Control for Best Scanning Linearity

Use our new digital control electronics with DDL (Dynamic Digital Linearization) to increase linearity and effective bandwidth in scanning applications by up to 1000-fold (see p. 6-16). By virtually eliminating tracking errors, DDL also increases the usable travel range.

### Working Principle and Reliability

P-540-series stages are equipped with the award-winning PICMA® piezo drives, integrated into a sophisticated, single-module, parallel-kinematics, flexure guiding system. The wire-EDM-cut flexures are FEA modeled for zero stiction, zero friction and exceptional guiding precision. The ceramic-encapsulated PICMA® drives are more robust than conventional piezo actuators, featuring superior lifetime and performance in both dynamic and static applications.

Because guidance, actuators and sensors are all frictionless and maintenance-free, these nanopositioning systems achieve outstanding levels of reliability.

### Notes

See the “Selection Guide” on p. 2-14 *ff.* for comparison with other nanopositioning systems.

### Ordering Information

**P-541.ZCD**  
Microscopy Z-Nanopositioning Stage, 100 µm, Direct Metrology, Capacitive Sensors

**P-541.TCD**  
Microscopy Z, Tip/Tilt Nanopositioning Stage, 100 µm / 1 mrad, Parallel Metrology, Capacitive Sensors

**Version with Strain Gauge Sensors see p. 2-48**



**P-541.ZSL**  
Microscopy Z Nanopositioning Stage, 100 µm, Strain Gauge Sensors

**P-541.TSL**  
Microscopy Z Tip / Tilt Nanopositioning Stage, 100 µm / 1 mrad, Strain Gauge Sensors

**Vacuum Versions Available.**

**P-540-Series XY Stages see p. 2-62**



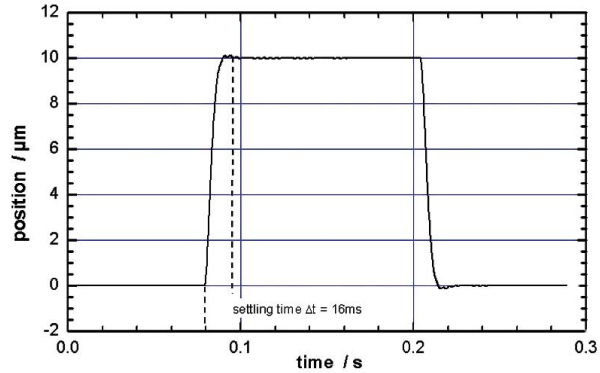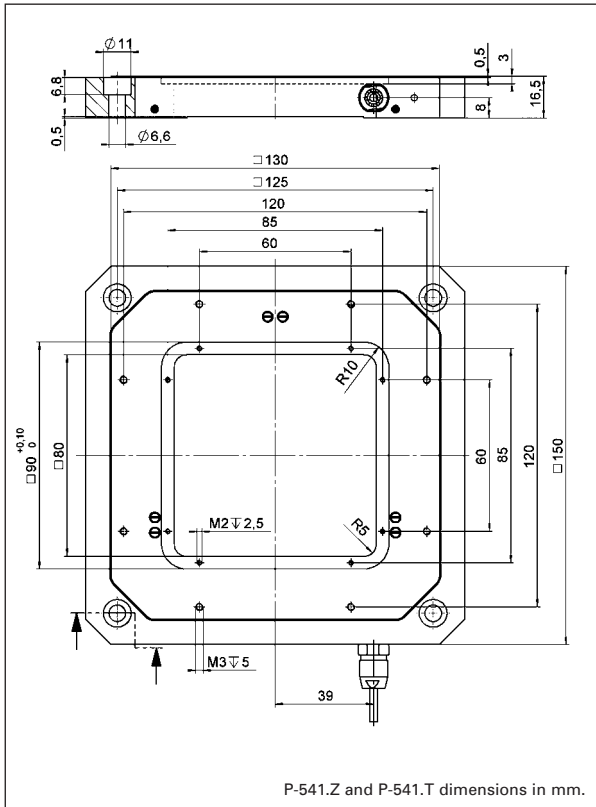

The settling time of a P-541.ZCD stage is approx. 16 ms for a 10  $\mu\text{m}$  step within a positioning accuracy of 1%.



For high-resolution Z-positioning systems with up to 460  $\mu\text{m}$  travel see the P-721 and P-725 PIFOC® nanofocusing drives (see p. 2-20 and p. 2-22).

## Technical Data

Models	P-541.ZCD	P-541.TCD	Units	Notes see p. 2-84
Active axes	Z	Z, $\theta_x$ , $\theta_y$		
Min. open-loop travel -20 to 120 V	150	150 (1.5 mrad)	$\mu\text{m}$	A2
Closed-loop travel	100	100 (1 mrad)	$\mu\text{m}$	A5
Integrated feedback sensor	capacitive	capacitive		B
* Closed-loop / open-loop resolution	0.8 / 0.2	0.8 / 0.2 (80 nrad / 20 nrad)	nm	C1
Closed-loop linearity (typ.)	0.03	0.03	%	
Repeatability		<5		nm
Push force capacity	100	100	N	D3
Pull force capacity	20	20	N	D3
Max. load	20	20	N	D4
Electrical capacitance	6.75	6.75	$\mu\text{F} \pm 20\%$	F1
** Dynamic Operating Current Coefficient	8.5	8.5 (Z)	$\mu\text{A}/(\text{Hz} \times \mu\text{m})$	F2
Resonant frequency unloaded	410	410 (Z)	$\text{Hz} \pm 20\%$	G2
Resonant frequency @ 185 g load	300	300 (Z)	$\text{Hz} \pm 20\%$	G3
Operating temperature	-20 to 80	-20 to 80	$^{\circ}\text{C}$	H2
Voltage & sensor connection	D	D		J1/J2
Body material	Al	Al		L
Recommended amplifier/controller (codes explained p. 2-17)	H, F, L	K		

\* For calibration information see p. 2-8. Resolution of PI piezo nanopositioning systems is not limited by friction or stiction. Noise equivalent motion with E-710, E-750, E-503 controllers / amplifiers.

\*\* Dynamic Operating Current Coefficient of linear axes is in  $\mu\text{A}$  per Hz and  $\mu\text{m}$ .

Example: Sinusoidal scan of 10  $\mu\text{m}$  at 10 Hz requires approximately 0.85 mA drive current.

# P-541.2CD · P-542.2CD

## Low-Profile, XY Piezo Scanning Microscopy Stages with Parallel Metrology



P-540-series nanopositioning stages feature a very low profile of 16.5 mm, a large 80 x 80 mm aperture and deliver highly accurate motion with sub-nanometer resolution.

- **Low Profile for Easy Integration: 16.5 mm**
- **Parallel Kinematics and Optional Parallel Metrology for Fast Response and Superior Multi-Axis Precision**
- **Choice of Strain Gauge (Lower Cost) or Capacitive Sensors (Higher Performance)**
- **To 200 x 200  $\mu\text{m}$  Travel Range**
- **Direct Drive Version for High-Speed Positioning & Scanning**
- **80 x 80 mm Clear Aperture**
- **PICMA® High-Performance Piezo Actuators for Superior Lifetime**

### Low Profile, Optimized for Microscopy Applications

P-541/P-542 nanopositioning and scanning stages are designed for easy integration into high-resolution microscopes. They feature a very low profile of 16.5 mm, a large 80 x 80 mm aperture, and offer highly accurate motion with sub-nanometer resolution.

#### Application Examples

- Optical trapping
- Scanning microscopy
- Mask & wafer alignment
- Scanning interferometry
- Surface metrology
- Biotechnology
- Micromanipulation

### Choice of Drives: Long Range or High-Speed Direct Drive

A variety of models are offered to suit a large range of applications: lever-amplified XY systems with 100 and 200  $\mu\text{m}$  travel and direct-driven XY scanners with 45  $\mu\text{m}$  travel. Their high resonant frequencies of 1.5 kHz in both axes allow for faster step response and higher scanning rates, needed for example in single-molecule microscopy, or in other time-critical applications. Z stages and Z-tip/tilt stages are also available (see p. 2-50).

### Higher Precision Through Parallel Kinematics/-Metrology with Capacitive Feedback Sensors

P-540-series XY piezo positioning stages feature a single-

module, parallel-kinematics design with all actuators operating on one central platform and no moving cables to cause microfriction. Advantages over serial kinematics setups are a lower profile, reduced inertia and better, axis-independent dynamics.

Capacitive sensors measure the actual distance between the fixed frame and the moving part of the stage directly and include any flex or other errors in the drive train—from the actuator through the lever and flexures to the platform—in the measurement. This results in higher motion linearity, long-term stability, phase fidelity, and—because external disturbances are seen by the sensor immediately—a stiffer, faster-responding servo-loop. See p. 2-4 *ff.* and p. 5-2 *ff.* for more information.

With parallel direct metrology, all capacitive sensors measure the position of the same moving platform against the same stationary reference (the fixed frame). This means that all motion is inside the servo-loop, no matter which actuator may have caused it, resulting in superior multi-axis precision. Advantages include: higher resolution in non-diffraction-limited imaging techniques (NSOM, etc.) and reduced blurring of edges in high-speed microscopy applications.

### Dynamic Digital Control for Best Scanning Linearity

Use our new digital control electronics with DDL (Dynamic Digital Linearization) to increase linearity and effective bandwidth in scanning applications by up to 1000-fold (see catalog E-710). By virtually eliminating tracking errors,

#### Ordering Information

**P-541.2DD**  
Microscopy XY Nanopositioning & Scanning Stage, High-Speed Direct Drive, 45 x 45  $\mu\text{m}$ , Parallel Metrology, Capacitive Sensors

**P-541.2CD**  
Microscopy XY Nanopositioning & Scanning Stage, 100 x 100  $\mu\text{m}$ , Parallel Metrology, Capacitive Sensors

**P-542.2CD**  
Microscopy XY Nanopositioning & Scanning Stage, 200 x 200  $\mu\text{m}$ , Parallel Metrology, Capacitive Sensors

**Versions with Strain Gauge Sensors see p. 2-60**

**P-541.2SL**  
Microscopy XY Nanopositioning & Scanning Stage, 100 x 100  $\mu\text{m}$ , Strain Gauge Sensors

**P-542.2SL**  
Microscopy XY Nanopositioning & Scanning Stage, 200 x 200  $\mu\text{m}$ , Strain Gauge Sensors


**Vacuum Versions Available.**

**P-540-Series Z-Tip/Tilt Stages see p. 2-50**

DDL also increases the usable travel range.

### Working Principle / Reliability

P-540-series stages are equipped with the award-winning PICMA® piezo drives, integrated into a sophisticated, single-module, parallel-kinematics, flexure guiding system. The wire-EDM-cut flexures are FEA modeled for zero stiction, zero friction and exceptional guiding precision. The ceramic-encapsulated PICMA® drives are more robust than conventional piezo actuators, featuring superior lifetime and performance in both dynamic and static applications. Because guidance, actuators and sensors are all frictionless and maintenance-free, these nanopositioning systems achieve outstanding levels of reliability.

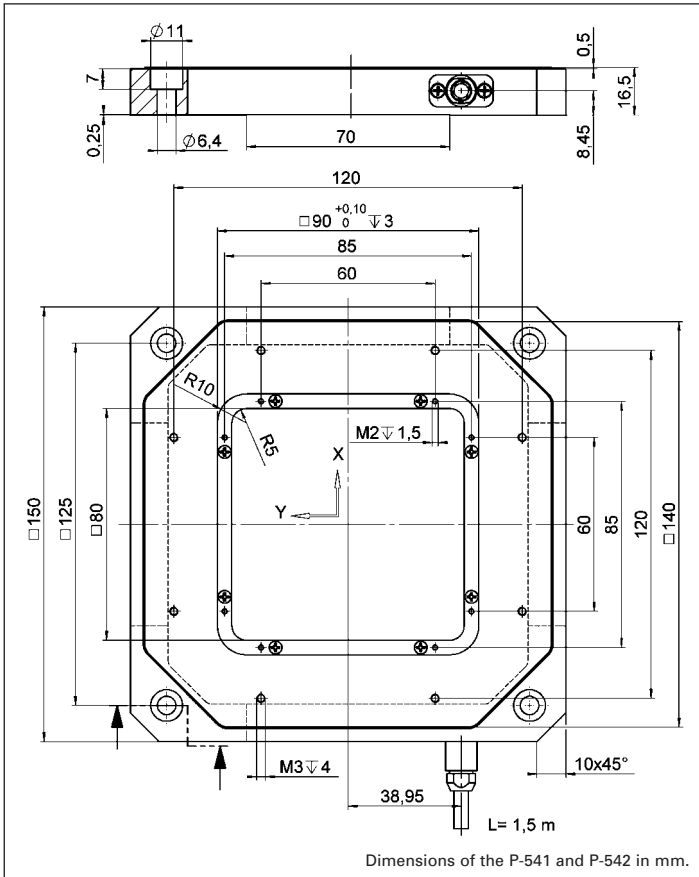

### Technical Data

Models	P-541.2CD	P-542.2CD	P-541.2DD	Units	Notes see p. 2-84
Active axes	XY	XY	XY		
Min. open-loop travel -20 to 120 V	150 x 150	250 x 250	60 x 60	μm	A2
Closed-loop travel	100 x 100	200 x 200	45 x 45	μm	A5
Integrated feedback sensor	capacitive	capacitive	capacitive		B
* Closed-loop / open-loop resolution	0.8 / 0.2	1.5 / 0.4	0.3 / 0.1	nm	C1
Closed-loop linearity (typ.)	0.03	0.03	0.03	%	
Repeatability	<5		<5	nm	
Push force capacity	100, 100	100, 100	200, 200	N	D3
Pull force capacity	20, 20	20, 20	20, 20	N	D3
Max. load	20	20	20	N	D4
Electrical capacitance (per axis)	6.75	7.5	28	μF ±20%	F1
** Dynamic Operating Current Coefficient (per axis)	8.5	4.8		μA/(Hz x μm)	F2
Resonant frequency unloaded	500, 500	370, 370	1500, 1500	Hz ±20%	G2
Resonant frequency @ 185 g load		250	1200, 1200	Hz ±20%	G3
Operating temperature	-20 to 80	-20 to 80	-20 to 80	°C	H2
Voltage & sensor connection	D	D	D		J1/J2
Body material	Al	Al	Al		L
Recommended amplifier/controller (codes explained p. 2-17)	H, F, L, K	H, F, L, K	H, F, L, K		

\* For calibration information see p. 2-8. Resolution of PI piezo nanopositioning systems is not limited by friction or stiction. Noise equivalent motion with E-710, E-750, E-503 controllers / amplifiers.

\*\* Dynamic Operating Current Coefficient in μA per Hz and μm. Example: Sinusoidal scan of 10 μm at 10 Hz with the P-542.2CD requires approximately 0.48 mA drive current.

# P-733.2DD · P-733.3DD

## Ultra-High-Speed, XY / XYZ Scanning Microscopy Stages with Parallel Metrology



P-733.3DD (left side) and P-733.2DD, high-speed, direct drive XY(Z) scanning stages are the fastest scanning stages with large aperture currently available (2.2 kHz resonant frequency!). Both units feature a footprint of only 100 x 100 mm. Compact disc for size comparison.

- Higher Speed Through Direct Drive
- Up to 2.2 kHz Resonant Frequency in X and Y
- 30 x 30 or 30 x 30 x 10  $\mu\text{m}$  Travel Range
- 100 Picometers Resolution
- Capacitive Sensors for Highest Linearity
- Parallel-Kinematics/Metrology for Enhanced Responsiveness / Multi-Axis Precision
- Active Runout Compensation
- 50 x 50 mm Clear Aperture
- Frictionless Precision Flexure Guiding System
- PICMA® High-Performance Piezo Drives

### Fastest Multi-Axis Systems w/Large Apertures

P-733.2DD / .3DD multi-axis piezo nanopositioning systems are the fastest ultra-high-precision, open-frame stages for

#### Application Examples

- Optical trapping
- Imaging (resolution enhancement)
- Scanning microscopy
- Surface structure analysis
- Biotechnology
- Atomic force microscopy
- Semiconductor test equipment
- Precision mask and wafer alignment
- Scanning interferometry
- Nanomanipulation
- Biophysics

microscopy. They provide a positioning and scanning range of 30 x 30 (x10)  $\mu\text{m}$  and are equipped with parallel-metrology capacitive position feedback for superior multi-axis linearity and repeatability.

The novel, high-stiffness, direct drive gives the systems resonant frequencies as high as 2.2 kHz (4 x that of other comparable systems), enabling millisecond scanning rates with sub-nanometer resolution.

### Low-Profile and Clear Aperture—Ideal for Microscopy

P-733 nanopositioning and scanning stages are designed for easy integration into high-resolution microscopes. They feature very low profiles, as low as 20 mm (0.8 inch), a 50 x 50 mm aperture, and offer highly accurate motion with sub-nanometer resolution.


### Higher Precision Through Parallel Kinematics/Metrology

P-733 piezo scanning stages feature a parallel-kinematics design with direct-measuring, non-contact capacitive position sensors (parallel, direct metrology).

These sensors make possible motion linearity to 0.03% with effective resolution in the sub-nanometer range. PI capacitive sensors are absolute-measuring, direct-metrology devices that boast very high bandwidth and exhibit no periodic errors.

Unlike conventional sensors, capacitive sensors measure the actual distance between the fixed frame and the moving part of the stage. They detect errors contributed by all components in the drive train—from the actuator through the flexures to the platform. This results in higher motion linearity, long-term stability, phase fidelity, and—because external disturbances are seen by the sensor immediately—a stiffer, faster-responding servo-loop. See p. 2-4 *ff.* and p. 5-2 *ff.* for more information.

Parallel kinematics means that all actuators act directly on the same moving platform leading to reduced size, inertia and the elimination of microfriction caused by moving cables. The advantages are enhanced dynamics, higher scanning rates, and better reproducibility.

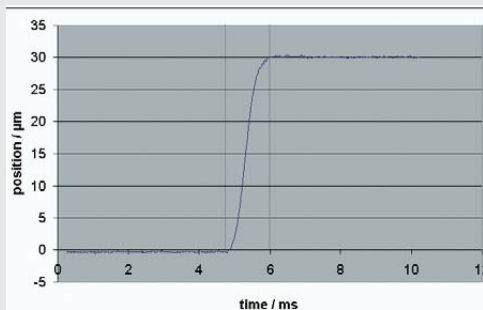

0 to 100 % travel in 1.36 msec, without overshoot. The P-733.2DD is the fastest, highest precision closed-loop XY scanning stage with large aperture currently available.

### Ordering Information

**P-733.2DD**  
XY High-Speed Direct Drive Piezo Scanning Stage, 30 x 30  $\mu\text{m}$ , Parallel Metrology, Sub-D Connector

**P-733.3DD**  
XYZ High-Speed Direct Drive Piezo Scanning Stage, 30 x 30 x 10  $\mu\text{m}$ , Parallel Metrology, Sub-D Connector

### Vacuum Versions Available


With parallel metrology, all sensors measure the position of the same moving platform against the same stationary reference (the fixed frame). This means that—in contrast to serial metrology—all motion is inside the servo-loop, no matter which actuator may have caused it, resulting in superior multi-axis precision (Active Trajectory Control).

### Dynamic Digital Control for Best Scanning Linearity

Use our new digital control electronics with DDL (Dynamic Digital Linearization) to increase linearity and effective bandwidth in scanning applications by up to 1000-fold (see p. 6-16). By virtually eliminating tracking errors, DDL also increases the usable travel range.

### Working Principle / Reliability

P-733 nanopositioning stages are equipped with the award winning PICMA® piezo drives,

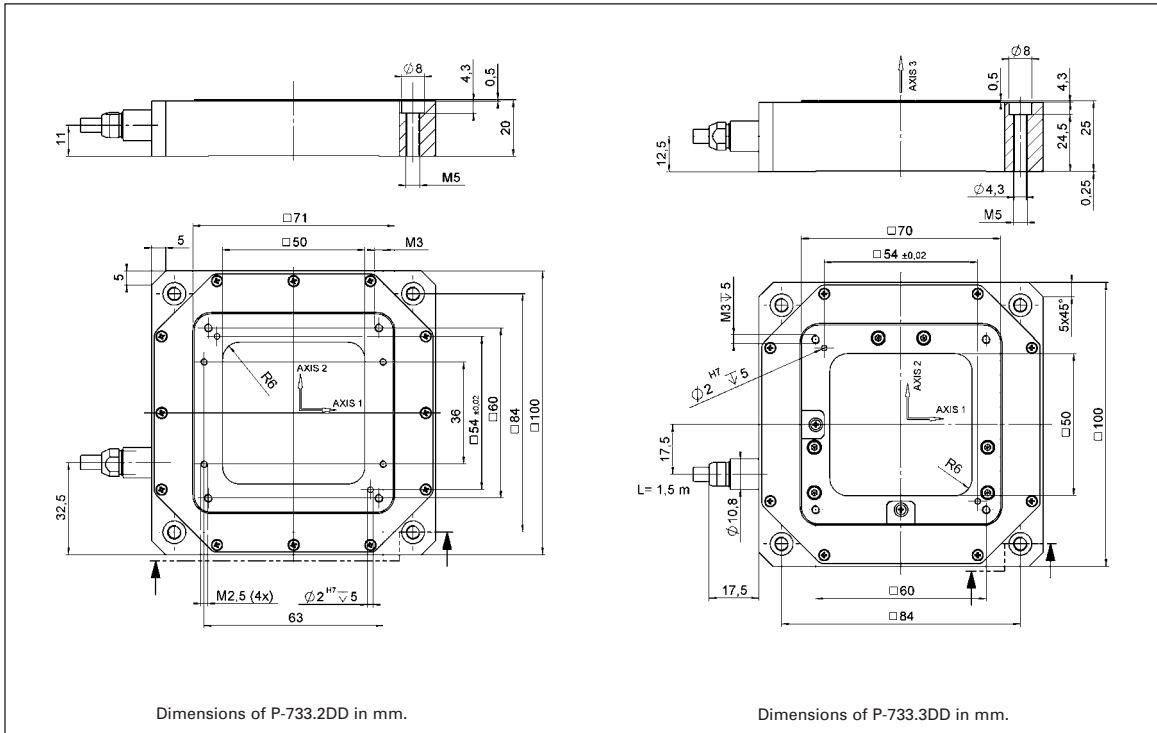

integrated into a sophisticated, single-module, parallel-kinematics, flexure guiding system. The wire-EDM-cut flexures are FEA modeled for zero stiction, zero friction and exceptional guiding precision. The ceramic-encapsulated PICMA® drives are more robust than conventional piezo actuators, featuring superior lifetime and performance in both dynamic and static applications.

Because guidance, actuators and sensors are all frictionless and maintenance-free, these nanopositioning systems achieve outstanding levels of reliability.

#### Vacuum Applications

PI offers versions specially designed for applications in ultra-high vacuum (see p. 2-66). A non-magnetizable version is available on request.

#### Notes

See the “Selection Guide” on p. 2-14 *ff.* for comparison with other nanopositioning systems.

#### Technical Data

Models	P-733.2DD	P-733.3DD	Units	Notes see p. 2-84
Active axes	X,Y	X,Y,Z		
Min. open-loop travel @ -20 to 120 V	33 x 33	33 x 33 x 14	$\mu\text{m} \pm 20\%$	A2
Closed-loop travel	30 x 30	30 x 30 x 10	$\mu\text{m}$	A5
Integrated feedback sensor	capacitive	capacitive		B
* Closed-loop / open-loop resolution	0.1 / 0.1	0.1 / 0.1	nm	C1
Closed-loop linearity (typ.)	0.03	0.03	%	
Full-range repeatability (typ.)	$\pm 1.0$	$\pm 1.0$	nm	C3
Stiffness (X,Y,Z)	20, 20, -	4, 4, 10	$\text{N}/\mu\text{m} \pm 20\%$	D1
Push / pull force capacity (in operating direction)	300 / 100	300 / 100	N	D3
Max. ( $\pm$ ) normal load	20	20	N	D4
Tilt $\theta_x, \theta_z$ (typ.)	3	3	$\mu\text{rad}$	E1
Electrical capacitance (X, Y, Z)	6.0, 6.0, -	6.0, 6.0, 4.4	$\mu\text{F} \pm 20\%$	F1
** Dynamic operating current (X, Y, Z) coefficient (DOCC)	25, 25,-	25, 25, 50	$\mu\text{A}/(\text{Hz} \times \mu\text{m})$	F2
Unloaded resonant frequency (X, Y, Z)	2230, 2230, -	1200, 1200, 1200	$\text{Hz} \pm 20\%$	G2
Resonant frequency @ 50 g load (X, Y, Z)	1800, 1800, -		$\text{Hz} \pm 20\%$	G3
Resonant frequency @ 200 g load (X, Y, Z)		530, 530, 635	$\text{Hz} \pm 20\%$	G3
Operating temperature range	-20 to 80	-20 to 80	$^{\circ}\text{C}$	H2
Voltage connection	sub-D, special	sub-D, special		J1
Sensor connection	sub-D, special	sub-D, special		J2
Weight (with cables)	525	635	$\text{g} \pm 5\%$	
Body material	Al	Al		L
Recommended amplifier/controller (codes explained p. 2-17)	H, F, L	H, F, L		

\* For calibration information see p. 2-8. Resolution of PI piezo nanopositioners is not limited by friction or stiction. The value given is noise equivalent motion with E-503 amplifier.

\*\* Dynamic Operating Current Coefficient in  $\mu\text{A}$  per Hz and  $\mu\text{m}$ . Example: Sinusoidal scan of 10  $\mu\text{m}$  at 100 Hz requires approximately 25 mA drive current.

# P-733

## High-Dynamics XY and XYZ Nanopositioning / Scanning Stages with Parallel Metrology



P-733.3CD XY scanning stage with active Z-axis  
P-733 stages are ideal for single-molecule microscopy

- Travel Ranges of 100 x 100 µm in X, Y and 10 µm in Z
- Direct Metrology with Capacitive Sensors
- Resolution to 0.1 nm
- Parallel Kinematics for Better Multi-Axis Accuracy and Dynamics
- Parallel Metrology for Active Trajectory Control
- 50 x 50 mm Clear Aperture for Transmitted-Light Applications

P-733 XY piezo driven stages are fast and highly accurate nanopositioning and scanning systems. They provide a positioning and scanning range of 100 x 100 (x10) µm together with sub-nanometer resolution. The 50 x 50 mm clear aperture is ideal for transmitted-light applications such as near-field scanning microscopy. The high-speed Z-axis of the P-733.3CD can actively compensate any out-of-plane, Z-axis deviation during XY motion.

### Application Examples

- Image processing / stabilization
- Scanning microscopy
- Metrology / Interferometry
- Biotechnology
- Semiconductor testing
- Mask / wafer positioning
- Micromanipulation
- Nanopositioning

### Capacitive Sensors for Higher Linearity

Capacitive sensors measure position directly and without physical contact. This makes them free of friction and hysteresis, a fact which, in combination with the positioning resolution of well under 1 nm, makes it possible to achieve very high levels of linearity. A further advantage of direct metrology with capacitive sensors is the high phase fidelity and the high bandwidth of up to 10 kHz. The closed-loop resolution is 0.3 nm for the X and Y axes and 0.2 nm for the optional Z-axis.

### Parallel Kinematics and Metrology for High Trajectory Fidelity

In a parallel kinematics multi-axis system, all actuators act directly on the same moving platform. This means that all axes see the same moved mass and can be designed with identical dynamic properties. Parallel kinematics systems have additional advantages

over serially stacked systems, including more-compact construction, no cumulative error from the different axes and no moving cables to cause friction that degrades resolution and linearity. Parallel kinematics systems can be operated with up to six degrees of freedom with low inertia and excellent dynamic performance.

Multiaxis nanopositioning systems equipped with both parallel kinematics and Direct Metrology are able to measure platform position in all degrees of freedom against a common, fixed reference: i.e. parallel metrology. In such systems, undesirable motion from one actuator in the direction of another (cross-talk) is detected immediately and actively compensated by the servo-loops. Known as Active Trajectory Control Concept, the parallel kinematics/parallel metrology combination can keep deviation from a trajectory to under a few nanometers, even in dynamic operation.

### Superior Lifetime

Highest possible reliability is assured by the use of award-winning PICMA® multilayer piezo actuators. PICMA® actuators are the only actuators on the market with ceramic-only insulation, which makes them resistant to ambient humidity and leakage-current failures. They are thus far superior to conventional actuators in reliability and lifetime.

### Wide Variety for Flexibility in Many Applications

For the highest dynamics applications, the P-733.2DD and P-733.3DD models with direct drive are available.

For Z-axis scanning applications, the P-733.ZCD version is

### Ordering Information

**P-733.2CD**  
XY Piezo Nanopositioning / Scanning Stage, 100 x 100 µm, Capacitive Sensors, Parallel Metrology, Sub-D Connector

**P-733.2CL**  
XY Piezo Nanopositioning / Scanning Stage, 100 x 100 µm, Capacitive Sensors, Parallel Metrology, LEMO Connector

**P-733.3CD**  
XYZ Piezo Nanopositioning / Scanning Stage, 100 x 100 x 10 µm, Capacitive Sensors, Parallel Metrology, Sub-D Connector

**P-733.3CL**  
XYZ Piezo Nanopositioning / Scanning Stage, 100 x 100 x 10 µm, Capacitive Sensors, Parallel Metrology, LEMO Connector

**P-733.2VL**  
XY Piezo Nanopositioning / Scanning Stage, 100 x 100 µm, Capacitive Sensors, Parallel Metrology, LEMO Connector, Vacuum Compatible down to 10<sup>-6</sup> hPa


**P-733.2VD**  
XY Piezo Nanopositioning / Scanning Stage, 100 x 100 µm, Capacitive Sensors, Parallel Metrology, Sub-D Connector, Vacuum Compatible down to 10<sup>-6</sup> hPa

**P-733.2UD**  
XY Piezo Nanopositioning / Scanning Stage, 100 x 100 µm, Capacitive Sensors, Parallel Metrology, Sub-D Connector, Vacuum Compatible down to 10<sup>-9</sup> hPa



Ask about custom designs.

available with a travel range of 100 µm.

For ultra-high-vacuum applications down to 10<sup>-9</sup> hPa, nanopositioning systems as well as comprehensive accessories, such as suitable feedthroughs, are available.

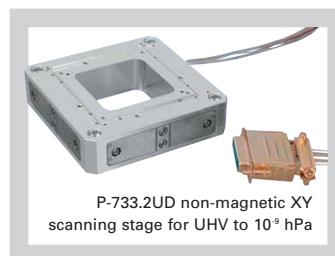

P-733.2UD non-magnetic XY scanning stage for UHV to 10<sup>-9</sup> hPa

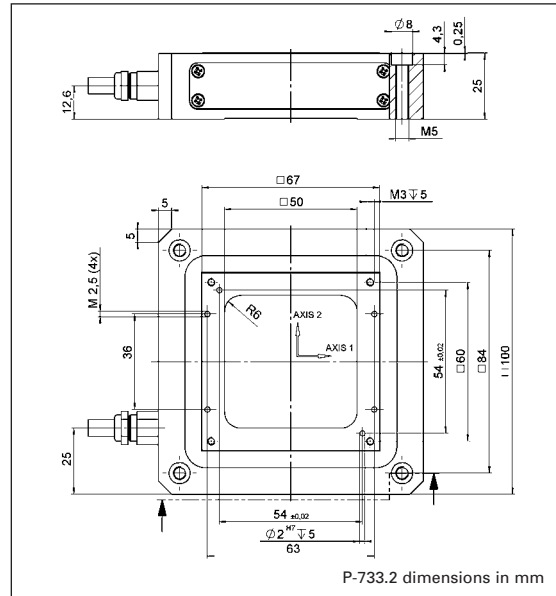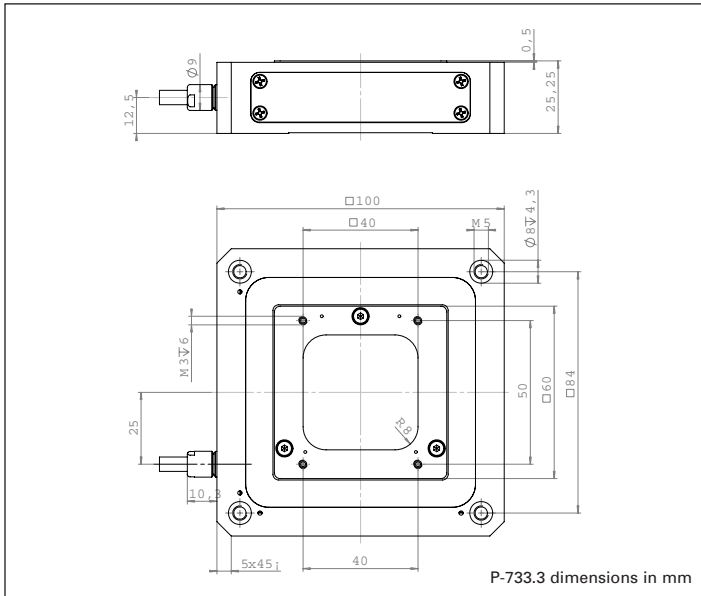

### Technical Data

	P-733.2CD	P-733.3CD	Units	Tolerance
Active axes	X, Y	X, Y, Z		
<b>Motion and positioning</b>				
Integrated sensor	Capacitive	Capacitive		
Open-loop travel, -20 to +120 V	115 x 115	115 x 115 x 12	μm	min.
Closed-loop travel	100 x 100	100 x 100 x 10	μm	
Open-loop resolution	0.2	0.2 in X, Y, 0.1 in Z	nm	typical
Closed-loop resolution	0.3	0.3 in X, Y, 0.2 in Z	nm	typical
Linearization	0.01	0.01 in X, Y, 0.03 in Z	%	typical
Repeatability	<2	<2 in X, Y, 1 in Z	nm	typical
Pitch	2	1	μrad	typical
Yaw	10	10	μrad	typical
<b>Mechanical properties</b>				
Stiffness in motion direction	1.4	1.4 in X, Y, 9.0 in Z	N/μm	±20 %
Unloaded resonant frequency	460	460 in X, Y, 1400 in Z	Hz	±20 %
Resonant frequency @ 120 g	340	340 in X, Y, 1060 in Z	Hz	±20 %
Resonant frequency @ 200 g	295	295 in X, Y, 610 in Z	Hz	
Push/pull force capacity in motion direction	300 / 100	300 / 100	N	Max.
Load	20	20	N	Max.
<b>Drive properties</b>				
Ceramic type	PICMA®	PICMA®		
Electrical capacitance	6.0 per axis	6.0 in X, Y, 3.3 in Z	μF	±20 %
Dynamic operating current coefficient	7.5 per axis	7.5 in X, Y, 40 in Z	μA/(Hz x μm)	±20 %
<b>Miscellaneous</b>				
Operating temperature range	-20 to +80	-20 to +80	°C	
Material	Al	Al		
Mass	580	675	g	±5 %
Cable length	1.5	1.5	m	±10 mm
Sensor connection	Sub-D Mix (CD-version) 4x LEMO (CL-version)	Sub-D Mix (CD-version) 6x LEMO (CL-version)		
Voltage connection	Sub-D Mix (CD-version) 2x LEMO (CL-version)	Sub-D Mix (CD-version) 3x LEMO (CL-version)		
Recommended controller/amplifier	H, F, L	H, F, L		

# P-733.Z

## High-Speed Piezo Z-Nanopositioning Stage with Direct Metrology



P-733.ZCD vertical piezo stage

- Travel Range 100  $\mu\text{m}$
- Fast: <8 ms Settling Time, 700 Hz Resonant Frequency
- Direct Metrology with Capacitive Sensors
- Closed-loop Resolution to 0.3 nm
- 50 x 50 mm Clear Aperture
- XY and XYZ Versions also Available
- Vacuum-Compatible Versions Available

P-733.Z piezo Z-stages offer a positioning and scanning range of 100  $\mu\text{m}$  with sub-nanometer resolution. The fast response and large clear aperture are ideal for applications such as fast focusing, scanning or 3D microscopy. The settling time of the P-733.Z is in the 10 ms range.

### Capacitive Sensors for Highest Accuracy

Capacitive sensors measure position directly and without physical contact. This makes them free of friction and hysteresis, a fact which, in combination with the positioning res-

olution of well under 1 nm, makes it possible to achieve very high levels of linearity. A further advantage of direct metrology with capacitive sensors is the high phase fidelity and the high bandwidth of up to 10 kHz. The resolution of the P-733.Z is better than 0.3 nm.

Capacitive position sensors measure the actual motion of the moving part relative to the stationary base (Direct Metrology). Errors in the drive train, actuator, lever arm or in guiding system do not influence the measurements. The result is exceptional motion linearity, higher long-term stability and a stiffer, more-responsive control loop, because external influences are immediately recognized by the sensor. The capacitive sensor non-linearity is typically less than 0.03%, the repeatability of the P-733.Z better than 2 nm.

### Superior Lifetime

Highest possible reliability is assured by the use of award-

winning PICMA® multilayer piezo actuators. PICMA® actuators are the only actuators on the market with ceramic-only insulation, which makes them resistant to ambient humidity and leakage-current failures. They are thus far superior to conventional actuators in reliability and lifetime.


### Large Variety of Models for a Broad Range of Applications

For scanning and positioning tasks in XY, the P-733.2CD and .3CD versions are available with a travel range of 100 x 100  $\mu\text{m}$ . For high-dynamics applications, the P-733.2DD and .3DD models can be offered with direct drive and reduced travel range.

For ultra-high-vacuum applications down to  $10^{-9}$  hPa, nanopositioning systems as well as comprehensive accessories, such as suitable feedthroughs, are available.

### Ordering Information

- P-733.ZCD**  
Z-Piezo Nanopositioning / Scanning Stage, 100  $\mu\text{m}$ , Capacitive Sensor, Sub-D Connector
- P-733.ZCL**  
Z-Piezo Nanopositioning / Scanning Stage, 100  $\mu\text{m}$ , Capacitive Sensor, LEMO Connector

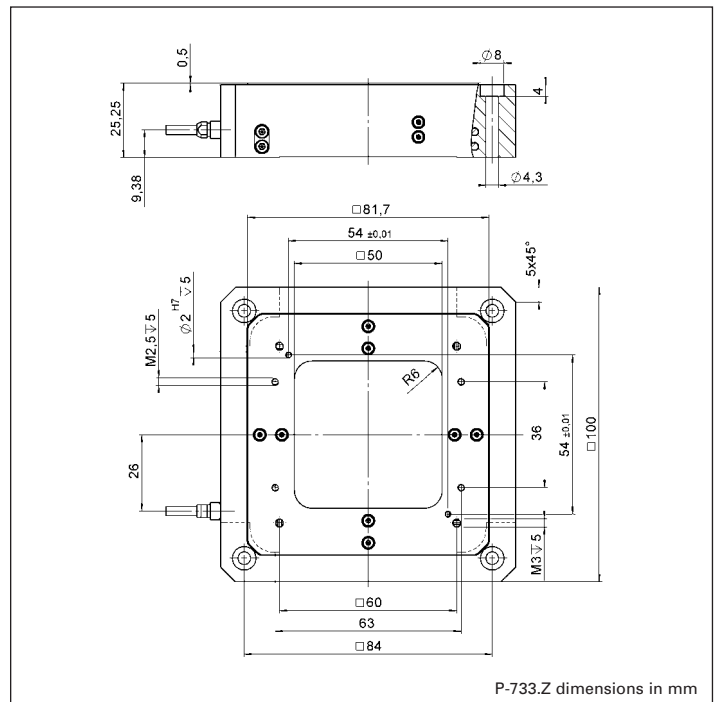


P-733.Z dimensions in mm

### Application Examples

- Scanning microscopy
- Confocal microscopy
- Mask/wafer positioning
- Surface inspection
- Nano-imprint
- Micromanipulation



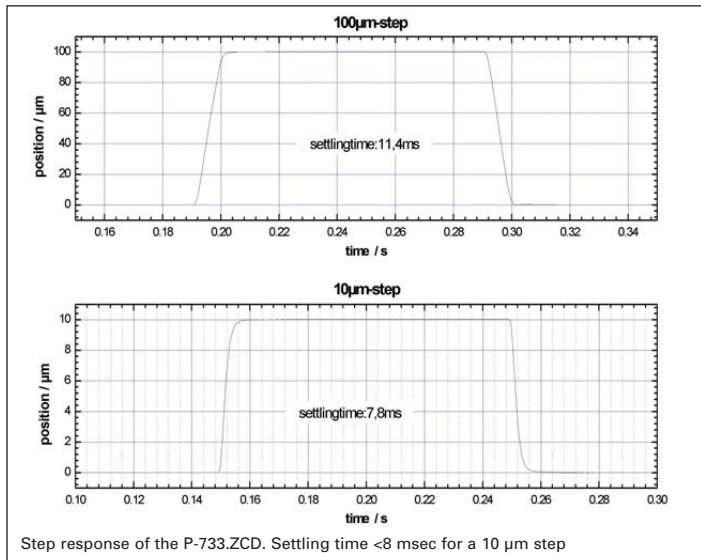


## Technical Data

	P-733.ZCD P-733.ZCL	Units	Tolerance
Active axes	Z		
<b>Motion and positioning</b>			
Integrated sensor	Capacitive		
Open-loop travel, -20 to +120 V	120	µm	min.
Closed-loop travel	100	µm	
Open-loop resolution	0.2	nm	typical
Closed-loop resolution	0.3	nm	typical
Linearity	0.03	%	typical
Repeatability	<2	nm	typical
Pitch	<5	µrad	typical
Yaw	<5	µrad	typical
<b>Mechanical properties</b>			
Stiffness in motion direction	2.5	N/µm	±20 %
Unloaded resonant frequency	700	Hz	±20 %
Resonant frequency @ 120 g	530	Hz	±20 %
Resonant frequency @ 200 g	415	Hz	±20 %
Load	20	N	Max.
<b>Drive properties</b>			
Ceramic type	PICMA®		
Electrical capacitance	6	µF	±20 %
Dynamic operating current coefficient	7.5	µA/(Hz x µm)	±20 %
<b>Miscellaneous</b>			
Operating temperature range	-20 bis 80	°C	
Material	Al		
Mass	580	g	±5%
Cable length	1.5	m	±10 mm
Sensor connection	Sub-D mix (CD-version) 2x LEMO (CL-version)		
Voltage connection	Sub-D mix (CD-version) 1x LEMO (CL-version)		
Recommended controller/amplifier	H, F, L		

# P-713 · P-714

## Low-Profile OEM XY Piezo-Scanners for Imaging Applications

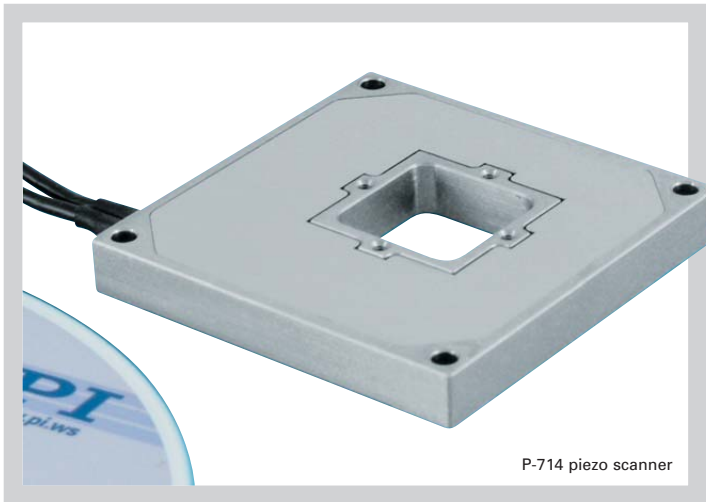

P-714 piezo scanner

- Ideal for Interlacing, CCD Resolution Enhancement
- Compact Size of only 45 x 45 x 6 mm
- Clear Aperture
- Highly Cost-Efficient Design
- 15 x 15 µm Travel Range
- Parallel-Kinematics Design for Higher Dynamics

P-713 / P-714 family piezo scanners and positioners with travel ranges of 15 x 15 µm feature especially compact designs. Ideal applications for the P-713 and P-714 are high-dynamics scanning or tracking tasks. Such tasks involve moving to specific positions in a small area (e.g. marked cells or CCD photosites) and from there following or performing motion with an amplitude of a few microns. The resonant frequency of up to over 2 kHz makes for settling times of a few millise-

conds, even after a full-range move, with closed-loop repeatability of under 5 nm.

### Flexibility

P-713 and P-714 nanopositioners are offered in different versions for different applications. The lowest-cost, basic version of the P-713 offers guiding accuracy in the motion plane of 50 µrad, a value generally good enough for interlacing tasks in scanning patterns of a few microns. For more demanding applications, the P-714 offers greater accuracy, typically 5 µrad or <10 nm absolute. If servo-control is required and no external position sensor is available, the P-714.2SL version is equipped with high-resolution strain gauge sensors (SGS) which offer nanometer-range resolution.

### Superior Lifetime

Reliability is assured by the use of award-winning PICMA® mul-

tilayer actuators, which are integrated into a sophisticated, single-module, flexure guiding system. The PICMA® actuators feature cofired ceramic encapsulation and thus offer better performance and reliability than conventional piezo actuators. The wire-EDM-cut flexures are FEA modeled to make them extremely precise. In addition they are maintenance-free and not subject to wear.

### Notes

See the "Piezo Drivers & Nanopositioning Controllers" section, p. 6-8 ff. for our comprehensive line of low-noise modular and OEM control electronics for computer and manual control.


See the "Selection Guide" on p. 2-14 ff. for comparison with other nanopositioning systems.

### Ordering Information

**P-713.20L**  
XY Piezo Scanner, 15 x 15 µm, SGS Sensor

**P-714.20L**  
XY Piezo Scanner, 15 x 15 µm, Improved Guiding Accuracy, Open-Loop



**P-714.2SL**  
XY Piezo Scanner, 15 x 15 µm, Improved Guiding Accuracy, SGS Sensors

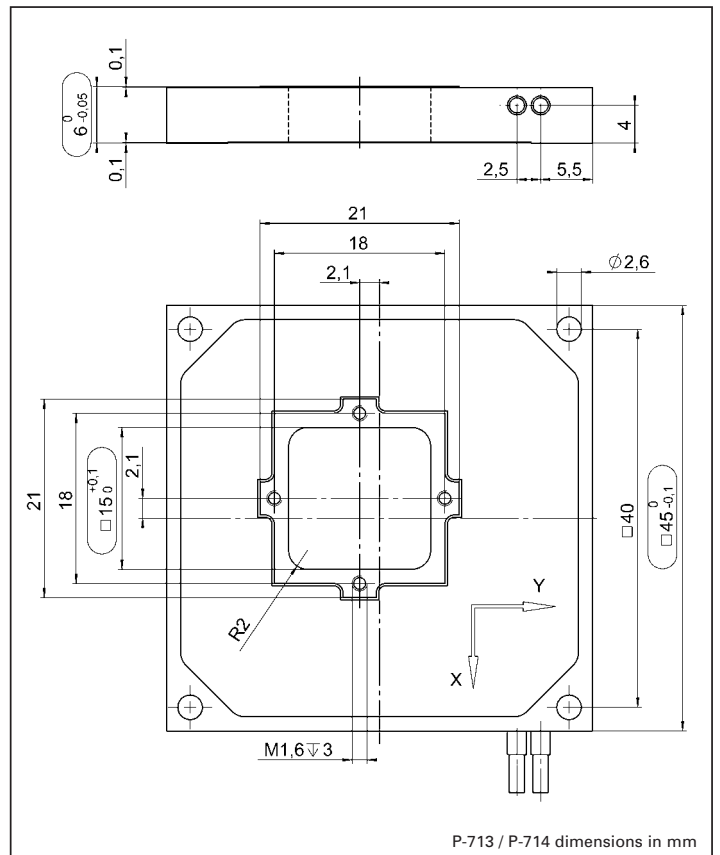

P-713 / P-714 dimensions in mm

### Application Examples

- Interlacing, image resolution enhancement
- Quality assurance
- Optical metrology
- Microscopy
- Imaging
- CCD camera technology

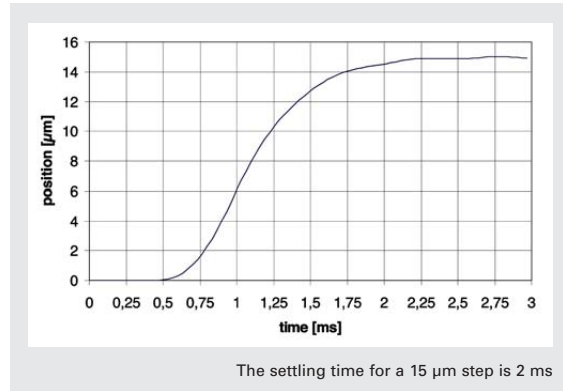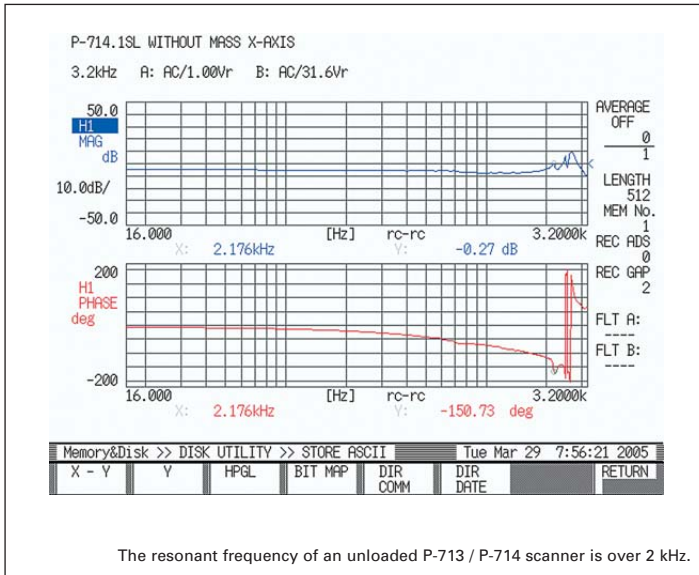

### Technical Data

Models	P-713.20L	P-714.20L	P-714.2SL	Units	Notes see p. 2-84
Active axes	X, Y	X, Y	X, Y		
Open-loop travel @ -20 to 120 V (X/Y)	20 / 20	20 / 20	25 / 20	µm ±20%	A2
Closed-loop travel	-	-	15 / axis	µm	A5
Integrated feedback sensor	-	-	SGS		B
* Open- / closed-loop resolution	0.1 / -	0.1 / -	1 / 0.1	nm	C1
Closed-loop linearity (typ.)	-	-	0.3	%	
Repeatability	-	-	<5	nm	
$\theta_x, \theta_y$	Typ. <1 µrad <5 µrad	Typ. <1 µrad <5 µrad	Typ. <1 µrad <5 µrad		
$\theta_z$	Typ. <30 µrad <50 µrad	Typ. <5 µrad <15 µrad	Typ. <5 µrad <15 µrad		
Stiffness	0.8	0.8	0.8	N/µm ±20%	D1
Max. normal load				kg	D4
Electrical capacitance	1.5 / axis	1.5 / axis	1.5 / axis	µF ±20%	F1
** Dynamic operating current coefficient (DOCC)	1.7 / axis	1.7 / axis	1.7 / axis	µA/(Hz x µm)	F2
Unloaded resonant frequency (X/Z)	2250	2250	2250	Hz ±20%	G2
Resonant frequency @ 20 g load	1310	1310	1310	Hz ±20%	G2
Resonant frequency @ 50 g load	1020	1020	1020		
Resonant frequency @ 100 g load	460	460	460	Hz ±20%	G2
Operating temperature range	-20 to 80	-20 to 80	-20 to 80	°C	H2
Voltage connection	VL	VL	VL		J1
Sensor connection	-	-	L		J2
Weight (with cable)	105	105	105	g ±5%	
Body material	S	S	S		L
Recommended amplifier/controller (codes explained p. 2-17)	A, G	A, G	D, H		

\* For calibration information see p. 2-8. Resolution of piezo nanopositioners is not limited by friction or stiction. Noise equivalent motion with E-503 amplifier.

\*\* Dynamic Operating Current Coefficient in µA per Hz and µm.

# P-363

## PicoCube® High-Speed, XY(Z) Piezo Stages for Nanotechnology, SPM, AFM



- Ultra-High-Performance Closed-Loop Scanner for AFM/SPM
- Compact Manipulation Tool for Nanotechnology
- Resonant Frequency 9.8 kHz
- Ultra-High-Precision Capacitive Feedback
- Parallel-Motion Metrology for Highest Linearity and Stability
- 50 Picometers Resolution
- 5 x 5 x 5 µm Travel Range
- Very Small Package
- Rugged Design

### Ultra-High-Performance Scanner/Manipulator for Nanotechnology

The P-363 PicoCube® XY/XYZ is an ultra-high-performance closed-loop piezo scanning system. Designed for AFM, SPM and nanomanipulation applications, it combines an ultra-low inertia, high-speed XY/XYZ piezo scanner with non-contact, direct-measuring, parallel-metrology capacitive

feedback capable of 50 picometers resolution. On top of being extremely precise, the PicoCube® system is also very small and rugged. Measuring only 30 x 30 x 40 mm (30 x 30 x 26 mm for XY version), it is easy to integrate in any scanning apparatus.

### AFM, STM, SPM, Nanolithography, Nanoimprinting, Nanometrology

The PicoCube® was specifically developed to overcome the limitations of the open-loop scanners currently available for STM, AFM and SPM. In addition to these applications, the PicoCube® is also the ideal scanning and manipulation tool for nanoimprinting, nanolithography, ultra-high-resolution, near-field, scanning optical microscopy and nano-surface-metrology applications.

### Higher Precision Through Parallel-Motion Metrology

The PicoCube® is based on a proprietary, ultra-fast, piezo-driven scanner design equipped with direct-measuring, capacitive position sensors (parallel metrology). Unlike conventional sensors, they measure the actual distance between the fixed frame and the moving part of the stage. This results in higher-motion linearity, long-term stability, phase fidelity, and—because external disturbances are seen by the sensor immediately—a stiffer, faster-responding servo-loop. See p. 2-4 *ff.* and p. 5-2 *ff.* for more information.

PI capacitive sensor electronics use the proprietary ILS (Integrated Linearization System) for enhanced linearity and are less sensitive to EMI than other high-resolution sensor systems.

With parallel metrology, all sensors measure the position of the same moving platform against the same stationary reference. This means that all motion is inside the servo-loop, resulting in superior multi-axis precision (Active Trajectory Control).

### Nanometer Accuracy in 1 Millisecond with 50-Picometer Resolution

PicoCube® systems provide resolution of 50 picometers and below. The ultra-fast XY/XYZ piezo drives offer resonant frequencies of 9.8 kHz in Z and >3 kHz in X and Y! The high resonant frequency and high-bandwidth capacitive feedback allow step and settle to 1% accuracy in as little as one millisecond.

### Rugged Design

In spite of its ability to move and position on an atomic scale, the PicoCube® boasts a

### Ordering Information

**P-363.3CD**  
PicoCube® XYZ Positioning- and Scanning System, 5 x 5 x 5 µm, Parallel-Motion Metrology, Capacitive Sensors, Sub-D Connectors

**P-363.3UD**  
Vacuum Version of P-363.3CD to 10<sup>-9</sup> hPa

**P-363.2CD**  
PicoCube® XY Positioning- and Scanning System, 5 x 5 µm, Parallel-Motion Metrology, Capacitive Sensors, Sub-D Connectors

**P-363.2UD**  
Vacuum Version of P-363.2CD to 10<sup>-9</sup> hPa

**P-363.3CL**  
PicoCube® XYZ Positioning- and Scanning System, 5 x 5 x 5 µm, Parallel-Motion Metrology, Capacitive Sensors, Lemo Connectors

**P-363.2CL**  
PicoCube® XY Positioning- and Scanning System, 5 x 5 µm, Parallel-Motion Metrology, Capacitive Sensors, Lemo Connectors

**Different travel ranges or dimensions on request.**


rugged design for real-world applications. For extra-high stability and reduced mass, the body is precision machined from heat-treated and stress-relieved titanium. The sophisticated frictionless design also ensures that the (moving) top plate protects the internal actuator/sensor unit from contamination.

### Control


The PicoCube® controller is based on the E-500, 19" rack-mount chassis with one E-509.CA3 servo-controller and three E-507.36 power amplifier modules. This controller is equipped with high-speed analog interfaces (0 to 10 V). An optional E-516 20-bit interface module (see page 6-26) is also available.

### Application Examples


- Atomic Force Microscopy
- Micromanipulation
- Biotechnology
- Nanomanipulation
- Nanoimprinting
- Nanometrology
- Nanolithography

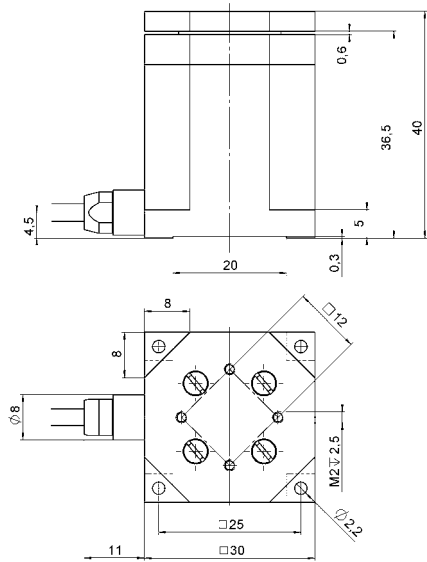


Dimensions of the P-363.3CD in mm. Removable top plate.

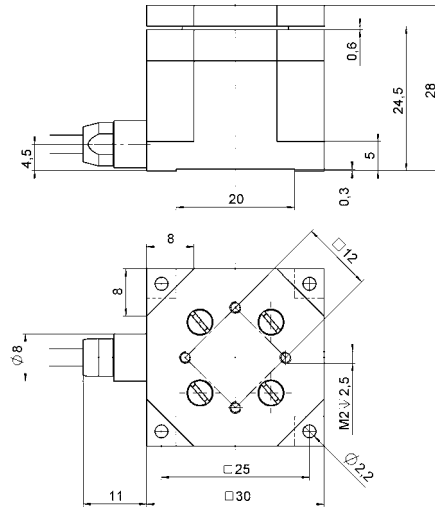

Dimensions of the P-363.2CD in mm. Removable top plate.

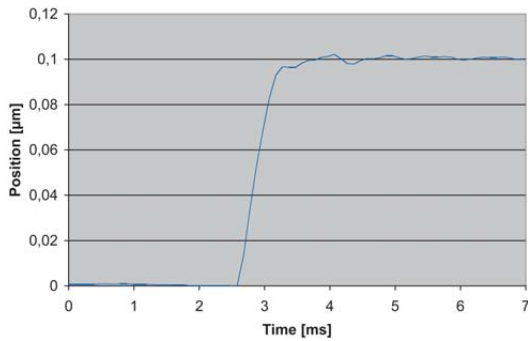

The P-363 settles to within 1 nm in 1 ms (100 nm step, X and Y motion; faster response in Z).

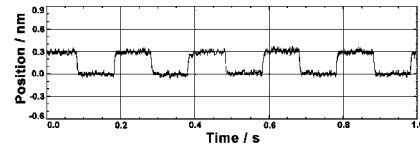

300 picometer steps (0.3 nm) performed with the P-363.

## Technical Data


Models	P-363.3CD	P-363.2CD	Units	Notes see p. 2-84
Active axes	X, Y, Z	X, Y		
Open-loop travel @ -250 to +250V	6, 6, 5.5	6, 6	µm	A5
Closed-loop travel	5, 5, 5	5, 5	µm	A2
Integrated feedback sensor	capacitive	capacitive		B
* Closed- / Open-loop resolution	0.1 / 0.05	0.1 / 0.05	nm	C1
Closed-loop linearity (X, Y, typ.)	1	1	nm	
Max. load	10	10	N	D4
Tilt, off-axis (typ.)	0.5	0.5	µrad	
Electrical capacitance	60 (X, Y); 110 (Z)	60	nF ±20%	F1
Unloaded resonant frequency (X, Y)	3.1	4.2	kHz ±20%	G2
Unloaded resonant frequency (Z)	9.8	–	kHz ±20%	G2
Resonant frequency @ 20 g load (X, Y)	1.5	2.1	kHz ±20%	G3
Operating temperature range	-40 to 120	-40 to 120	°C	H2
** Voltage and sensor connection	Sub-D special	Sub-D special		J1
Weight (with cables)	225	190	g ±5%	
Body material	Titanium	Titanium		L
Recommended amplifier/controller	E-500 rack, E-509, E-507.36 amplifier	E-500 rack, E-509, E-507.36 amplifier		

\* For calibration information see p. 2-8.  
Resolution of PI piezo nano-positioners is not limited by friction or stiction. The value given is noise equivalent motion with E-509.C3A controller and E-507.36 amplifier.  
\*\* P-363.xCL versions with Lemo connectors

## News & Applications

# World's First Digital Piezo Nanopositioning Controller on PCI Board

New E-761 digital piezo controller in the PCI-board format

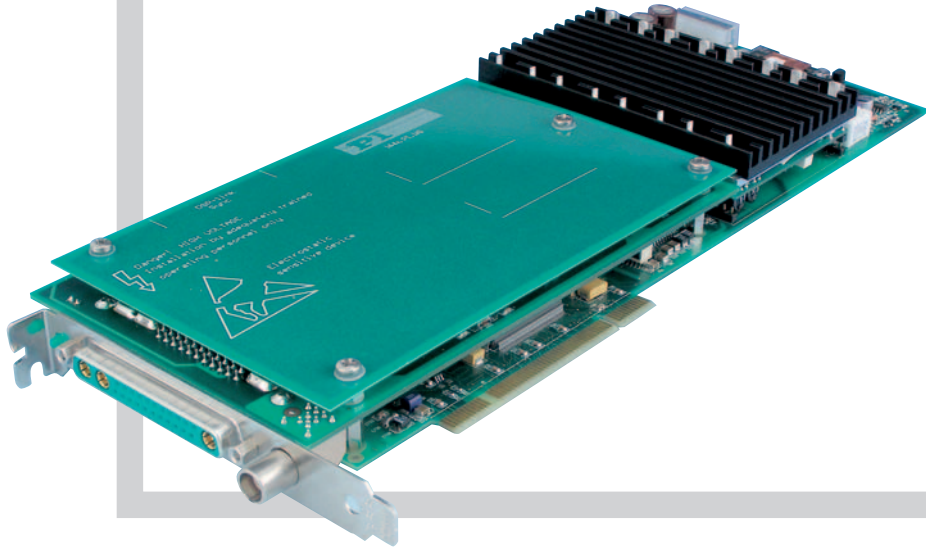

The PCI bus allows for very rapid communication and easy integration with devices such as frame grabbers – a feature which is very advantageous in real-time data acquisition applications or when operating multiple axes or controllers simultaneously.

The internal coordinate transformation means it is no problem to operate complex, multi-axis stages, such as three-axis X-Y- $\Theta_z$  stages or Z-tip-tilt platforms. In systems with parallel metrology, the E-761 can also automatically compensate undesired off-axis motion (active trajectory control), making it possible to attain motion accuracies in the sub-nanometer range.

Of course, this digital controller has the technical refinements you have come to expect from PI, such as 32-bit digital filters, 24-bit DAC resolution, multi-stage sensor and electronics linearization, plug & play ID-chip support and an extensive software support package.


### E-761 3-channel digital nanopositioning controller: very advanced control technology in cost-effective PCI format to control up to three logical axes

The E-761 digital nanopositioning controller offers very advanced control technology in the cost-effective PCI-board format. It is able to control nanopositioning systems with up to three logical axes and four piezo actuators.

Many of today's high-tech applications, such as imaging, metrology, scanning microscopy and surface analysis require a combination of high-speed motion control and high-resolution vision. They also require extremely fast data acquisition and precise synchronization between the imaging and motion control devices. Peripheral components with PCI bus interface are ideal for these tasks, because the PCI bus was designed to give high-bandwidth access to the microprocessor in the PC. With the E-761, PI intro-


duces the world's first fully digital piezoelectric nanopositioning controller on a PCI board.

P-541.2CD low-profile XY piezo stage with capacitive sensors: a perfect fit for the E-761

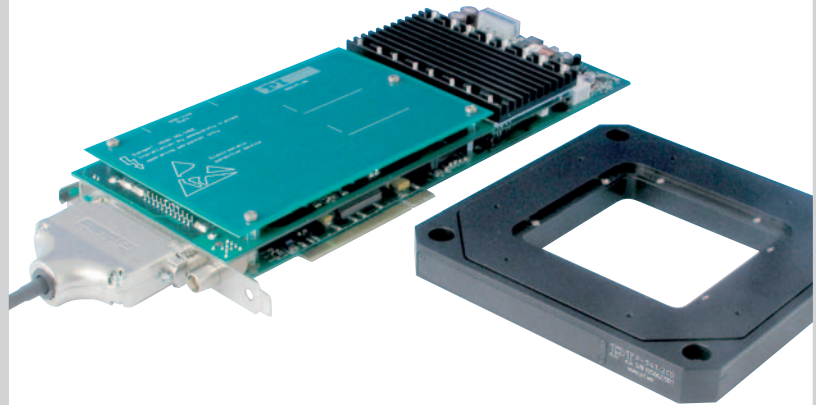

## News & Applications

# PIFOC® Revisited: P-737 Piezo-Z Microscopy Specimen Scanning Stage

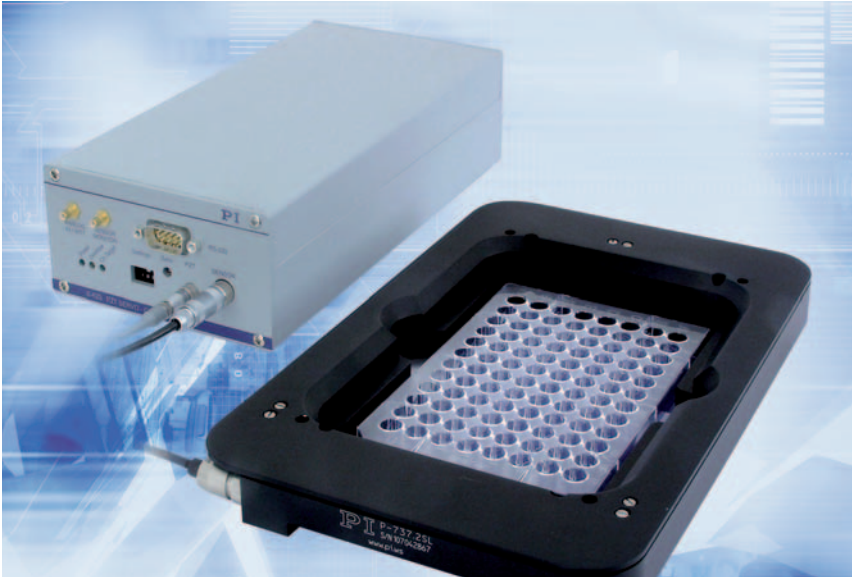


PIFOC® P-737 high-speed vertical positioning systems are designed for integration into motorized XY microscopy stages of leading manufacturers such as Märzhäuser. While the XY stage positions the sample, the piezo flexure-based P-737 moves the sample along the optical axis to adjust the focus quickly and precisely. Vertical stepping with nanometer precision takes only a few milliseconds. The large aperture is designed to accommodate a variety of specimen holders including slides or multiwell plates.

### High-speed Z-steps for fast focus control and imaging

The immediate response of the solid-state piezo drives enables rapid Z-steps with typically 10-to-20-times-faster step-and-settle than classical stepper motor drives. This leads to significantly higher image acquisition speed and throughput. The P-737 is available with travel ranges of 100 or 250  $\mu\text{m}$  and offers a choice of feedback sensors between strain gauge or non-contact, capacitive position sensors, depending on the accuracy requirements.

PI is the inventor of high-speed Piezo-Z objective steppers. Today the term PIFOC® is practically synonymous with all high-resolution vertical drives for microscope objectives. Depending on the application, it can be advantageous to adjust the sample instead of the objective. While no sample stage can beat the speed of the fastest PIFOC® objective positioners (due to the stiffer and more compact design), the P-737 stage allows for convenient Z-stack imaging with multiple objectives and very high throughput. In both cases the effect remains the same: the focal plane moves through the sample. This is why the P-737 is also marketed under the PIFOC® trademark.

Settling time is in the 10 ms range – here, 1  $\mu\text{m}$  steps

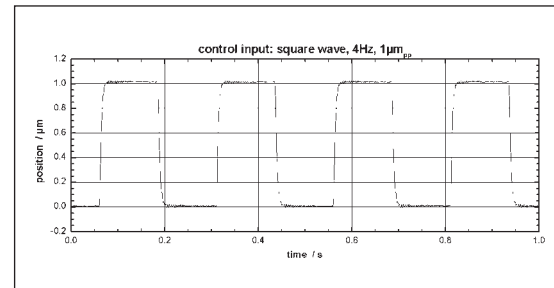

### Analog or Digital Systems

The P-737 together with the E-625 piezo-controller, offers a cost-effective system for high-speed, high-resolution positioning of microscopy specimens. The controller features a choice of a 20-bit digital interface or a broadband analog interface for the target position.

The capacitive-sensor-equipped version can also be operated with the sophisticated E-753 digital servo-controller. The main advantage of this controller is its adaptability to changing load conditions on the piezo stage. No matter what load is applied, the system can always provide an optimum combination of speed, settling time and resolution.



PIFOC® objective positioners and scanners are available with travel ranges of up to 400  $\mu\text{m}$ . QuickLock thread adapters enable fast mounting of the PIFOC® on the microscope and flexible replacement of objective

## News & Applications

New Multi-Axis Piezo Controller Provides Higher Precision, Flexibility

# Digitally Controlling the Nano-World – in up to 6 Axes

The new, E-712 multi-axis piezo-controller picks up where the successful E-710 left off. It features a faster processor, a real-time operating system and significantly higher servo update and sensor sampling rates to provide extremely precise coordinated motion in up to 6 degrees of freedom with nano-meter precision.



### Additional advantages:

- Versatile interfaces: TCP/IP for remote control through the internet; additional USB 2.0 and the classic RS-232. Available high-resolution analog inputs for direct control with high-bandwidth analog signals.
- 20-bit D/A converters make possible sub-nanometer position resolution, even over long travel ranges of  $>1000 \mu\text{m}$ , like those of the new P-629 PIHera® piezo nanopositioning stages.
- Available high-power amplifier modules, to supplement the integrated amplifiers for high-frequency scanning / tracking applications.
- Software – tunable, proportional-integral, digital servo-control with 2 notch filters allows operation of the piezomechanics closer to its resonant frequency.
- 600 MHz processor and sensor and servo update rates of up to 50 kHz assure faster updating of the position and control data – especially important for high-dynamics periodic motion. In addition, optional Digital Dynamic Linearization, ensures that tracking errors are reduced to the nanometer level even in high-dynamics scanning applications.
- Internal coordinate transformation for parallel-kinematics systems with user-friendly position commands in Cartesian coordinates.

The E-712 controller is available in 3- and 6-channel configurations.

Mature operating software with intuitive interfaces make operating the system almost child's play. No programming knowledge is needed, whether commissioning, optimizing system parameters, creating new user-defined motion profiles for the integrated wave generator, or recording data generated during the motion. The comprehensive package of supporting software, including LabVIEW drivers and DLLs, ensures easy integration in a variety of system environments.



## News & Applications

# PicoCube®: Reference-Class System for Nanotechnology



**New E-536 controllers enable even higher resolution and bandwidth with the minute P-363 PicoCube® piezo stages.**

The minute P-363 PicoCube®, together with its reference-class, ultra-low-noise E-536 driver/controller, are ideally suited for nanotechnology applications. They provide significantly higher dynamics, resolution and positional stability than previous multi-axis scanning stages. Nano-imprint lithography, scanning microscopy and biotechnology benefit from the extremely high resolution of up to 25 picometers over travel ranges of 5  $\mu\text{m}$  per axis.

### **Controllers: Optimized for Highest Resolution / Bandwidth**

Two controller versions are available. For high-speed scanning applications the E-536.3C high-power models featuring 100 watts per chan-

nel are recommended. The E-536.3CH ultra-low noise versions are optimized for highest positioning accuracy and resolution in the picometer range.

### **Compact 2-Axis and 3-Axis Stages**

The compact PicoCube® is available as XY and XYZ system. It is based on exceptionally robust, high-stiffness shear piezo drives and employs non-contact, direct-measuring, parallelmetrology capacitive sensors for position feedback.

The low-inertia drives allow for a resonant frequency of 10 kHz, important for high speed scanning applications. Measuring only 30 x 30 x 26 mm (XY version), it is easy to integrate in any scanning apparatus.

**For more information, see p. 28**

# Overview: Piezo Controllers and Drivers

## For Closed-Loop & Open-Loop Piezo Stages / Actuators, Digital & Analog Control

### E-710 Multi-Channel Digital Piezo Controllers with Dynamic Linearization



- For Piezo Nanopositioners with Capacitive Feedback Sensors
  - 3-, 4- & 6-Channel Versions, 32-Bit Digital Filters
  - Hardware and Firmware Linearization
  - Coordinate Transformation for Parallel Kinematics /Parallel Metrology Systems
  - Optional Dynamic Digital Linearization (Firmware Option) Improves Scanning Linearity
  - ID Chip Support for Automatic Calibration
  - Optional High-Speed Parallel I/O or DSP Interface
  - Optional Analog Input
  - GCS (General Command Set) Compatible
- See Online Datasheet Available at [www.pi.ws](http://www.pi.ws)**


### E-712 Modular Digital Piezo Controller System, up to 6 Axes



- Digital Controller of the Newest Generation: 600 MHz Processor; up to 50 kHz Sensor Update Rate; Highly Stable 20-bit D/A Converters
- Real-Time Operating System for Excellent Trajectory Control
- Modular Design for Greatest Flexibility in Meeting Custom Requirements
- Auto-Calibrating from Stage ID-Chip for Interchangeability of Controller and Mechanics w/o Recalibration
- Flexible Interfacing: USB, RS-232 & Ethernet
- Optional High-Bandwidth Analog Inputs and Outputs

**See p 32 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)**

### E-753 Digital High-Speed Piezo Controller for Single-Channel Nanopositioning Systems

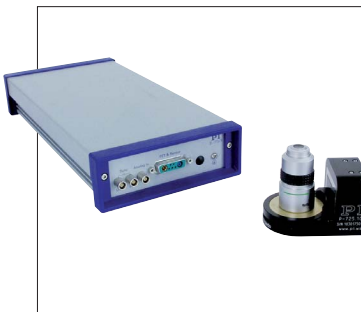

- Next Generation Digital Controller Provides Higher Flexibility, Accuracy and Speed
- 100 kHz Sensor Sampling; 32-bit Floating Point DSP; 24-bit Low-Noise D/A Converters
- Automatic Calibration Data Coaching from Stage ID-Chip
- Additional High-Bandwidth Analog Control Input / Sensor Input
- Digital I/O Lines for Triggering Tasks
- Extensive Software Support
- For Nanopositioning Systems with Capacitive Sensors

**See Online Datasheet Available at [www.pi.ws](http://www.pi.ws)**

### E-500 – E-501 Modular Piezo Control Systems (HVPZT & LVPZT)




- For High- and Low-Voltage Piezo Systems
- 19- and 9.5-Inch Chassis
- Optional Position Servo-Control Modules
- Optional Display and Digital Interface Module

**See Online Datasheet Available at [www.pi.ws](http://www.pi.ws)**

# Piezo Controllers and Drivers

## E-610, E-625, E-831 Piezo Amplifiers / Controller Modules with Optional Computer Interface

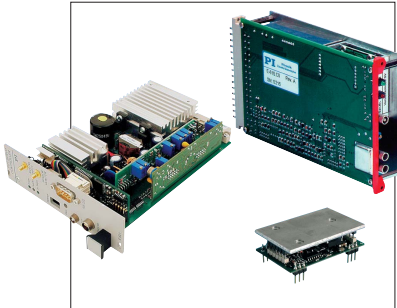

- Analog & Optional 20-Bit High Speed Digital Interface
- Single Channel, Network Capability with up to 12 Channels
- For Capacitive, LVDT and Strain Gauge Position Sensors
- Open-Loop and Closed-Loop Versions
- Stand Alone or Rack-Mount Operation

See Online Datasheet Available at [www.pi.ws](http://www.pi.ws)


## E-755 Digital Controller for NEXLINE® Nanopositioning Linear Drives



- Special Control Algorithms for NEXLINE® Nanopositioning Linear Drives
- 32-Bit Digital Filters
- 24-Bit DAC Resolution
- Fully Programmable Low-Pass and Notch Filters
- Non-Volatile User Settings and Last-Position Data
- PI General Command Set Compatible

See Online Datasheet Available at [www.pi.ws](http://www.pi.ws)


## E-625 and E-665 Piezo Amplifiers & Servo-Controllers

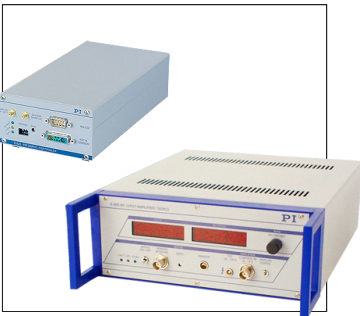

- Analog & 20-Bit High-Speed Digital Interface
- SGS, LVDT and Capacitive Sensor Servo-Control
- Wave Table for User-Defined Curves
- Network Capability with up to 4 or up to 12 Channels
- GCS (General Command Set) Compatible

See Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

## E-761 Digital Piezo Controller PCI Board for Piezo Stages

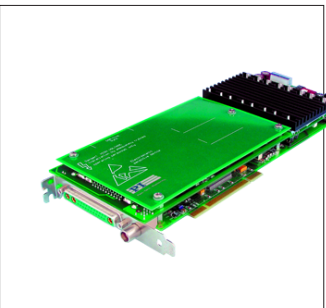

- For Piezo Stages with Capacitive Sensors
- High-Speed PCI Interface
- 3 Logical Axes, 4 Piezo Amplifiers
- Additional High-Bandwidth Analog Interface
- 32-Bit Digital Filters
- 24-Bit Ultra-Low-Noise DAC Converters
- Coordinate Transformation for Parallel-Kinematics / Parallel-Metrology Systems

See p. 30 or Online Datasheet Available at [www.pi.ws](http://www.pi.ws)

### Program Overview

- Piezo Ceramic Actuators & Motors
- Piezo Nanopositioning Systems and Scanners
- Active Optics / Tip-Tilt Platforms
- Capacitive Nanometrology Sensors
- Piezo Electronics: Amplifiers and Controllers
- Hexapod 6-Axis Positioners / Robots
- Micropositioning Stages & Actuators
- Photonics Alignment Systems, Solutions for Telecommunications
- Motor Controllers
- Ultrasonic Linear Motors

Request or download the complete  
PI Nanopositioning & Piezo Actuator  
Catalog



Call or go to: <http://www.pi.ws>



## Headquarters

### GERMANY

#### Physik Instrumente (PI) GmbH & Co. KG

Auf der Römerstr. 1  
D-76228 Karlsruhe/Palmbach  
Tel: +49 (721) 4846-0  
Fax: +49 (721) 4846-100  
info@pi.ws · www.pi.ws

#### PI Ceramic GmbH

Lindenstr.  
D-07589 Lederhose  
Tel: +49 (36604) 882-0  
Fax: +49 (36604) 882-25  
info@piceramic.de  
www.piceramic.de

## Subsidiaries

### USA (East) & CANADA

#### PI (Physik Instrumente) L.P.

16 Albert St.  
Auburn, MA 01501  
Tel: +1 (508) 832 3456  
Fax: +1 (508) 832 0506  
info@pi-usa.us  
www.pi-usa.us

### USA (West) & MEXICO

#### PI (Physik Instrumente) L.P.

5420 Trabuco Rd., Suite 100  
Irvine, CA 92620  
Tel: +1 (949) 679 9191  
Fax: +1 (949) 679 9292  
info@pi-usa.us  
www.pi-usa.us

### JAPAN

#### PI Japan Co., Ltd.

Akebono-cho 2-38-5  
Tachikawa-shi  
J-Tokyo 190  
Tel: +81 (42) 526 7300  
Fax: +81 (42) 526 7301  
info@pi-japan.jp  
www.pi-japan.jp

#### PI Japan Co., Ltd.

Hanahara Dai-ni Building, #703  
4-11-27 Nishinakajima,  
Yodogawa-ku, Osaka-shi  
J-Osaka 532  
Tel: +81 (6) 6304 5605  
Fax: +81 (6) 6304 5606  
info@pi-japan.jp  
www.pi-japan.jp

### CHINA

#### Physik Instrumente (PI Shanghai) Co., Ltd.

Building No. 7-301  
Longdong Avenue 3000  
201203 Shanghai, China  
Tel: +86 (21) 687 900 08  
Fax: +86 (21) 687 900 98  
info@pi-china.cn  
www.pi-china.cn

### GREAT BRITAIN

#### PI (Physik Instrumente) Ltd.

Lambda House  
Batford Mill  
GB-Harpenden, Hertfordshire  
AL5 5BZ  
Tel: +44 (1582) 764 334  
Fax: +44 (1582) 712 084  
info@physikinstrumente.co.uk  
www.physikinstrumente.co.uk

### FRANCE

#### PI France S.A.S.

32 rue Delizy  
F-93694 Pantin Cedex  
Tel: +33 (1) 481 039 38  
Fax: +33 (1) 481 009 66  
info@pi-france.fr  
www.pi-france.fr

### ITALY

#### Physik Instrumente (PI) S.r.l.

Via G. Marconi, 28  
I-20091 Bresso (MI)  
Tel: +39 (02) 665 011 01  
Fax: +39 (02) 665 014 56  
info@pionline.it  
www.pionline.it