


# Microscopy and Industrial Applications

## What is your Application?

### 2 (3) Photon Fluorescence Microscopy

In 2 or 3-photon fluorescence microscopy, two photons excite visible fluorescence in deep tissue layers in the infrared range. By using infrared emissions, a deeper penetration is achieved, this is why a respectively large travel range of the focus drive is required. For example, cell organs, metabolic processes, etc. were examined. The users generate Z image stacks.

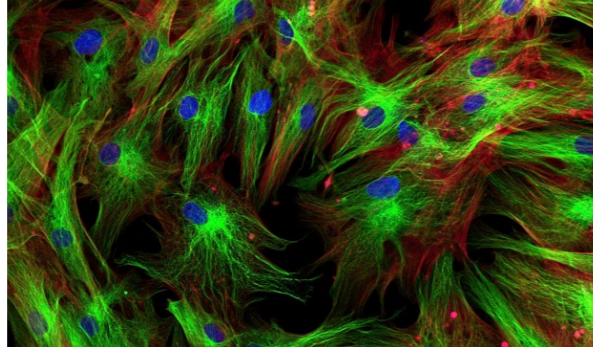


### Digital Slide Scanning

Today most microscopy routine diagnostics on histological slides (histology, pathology, hematology, gynecology) is done by automated screening devices. These devices scan the sample slides in a structured way at a very high speed. To adjust the z-focus and to accommodate to z-height differences of the sample during the XY-scanning the Z-focus drive has to quickly and precisely adjust in travel ranges typically between several hundred micrometer to several millimeter, depending on the device design.

### Wafer Inspection

As semiconductor devices move towards higher integration and density with advances in semiconductor materials and processes, the critical defect size will decrease. There are several factors that affect the topology of the wafer surface: The surface tension induced by IC manufacturing and the surface roughness. To quickly obtain a high-quality surface image, defect inspection tools must operate at high resolutions and with high dynamics. Precise and fast autofocus ability is critical.

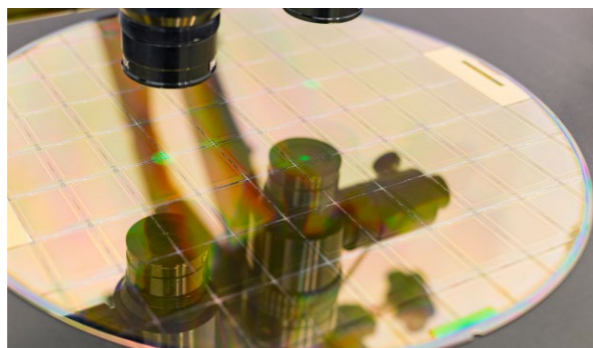

## Headquarters

### GERMANY

**Physik Instrumente (PI) GmbH & Co. KG**  
Auf der Roemerstrasse 1  
76228 Karlsruhe  
Phone +49 721 4846-0  
Fax +49 721 4846-1019  
info@pi.ws  
www.pi.ws

### PI miCos GmbH

Freiburger Strasse 30  
79427 Eschbach  
Phone +49 7634 5057-0  
Fax +49 7634 5057-99  
info@pimicos.com  
www.pi.ws

### PI Ceramic GmbH

Lindenstrasse  
07589 Lederhose  
Phone +49 36604 882-0  
Fax +49 36604 882-4109  
info@piceramic.com  
www.piceramic.com

## ACS Motion Control

### ISRAEL

**ACS Motion Control Ltd.**  
Ramat Gabriel Industrial Park  
1 Hataasia St.  
Migdal HaEmek, 2307037  
POB 984  
Phone +972-4-6546440  
Fax +972-4-6546443  
info@acsmotioncontrol.com  
www.acsmotioncontrol.com



© Physik Instrumente (PI) GmbH & Co. KG

All contents, including texts, graphics, data etc., as well as their layout, are subject to copyright and other protective laws. Any copying, modification or redistribution in whole or in parts is subject to a written permission of PI.

Although the information in this document has been compiled with the greatest care, errors cannot be ruled out completely. Therefore, we cannot guarantee for the information being complete, correct and up to date. Illustrations may differ from the original and are not binding. PI reserves the right to supplement or change the information provided without prior notice.

Follow us on:



## PI Subsidiaries

### USA (East) & CANADA

**PI (Physik Instrumente) L.P.**  
Auburn, MA 01501  
www.pi-usa.us

### USA (San Francisco Bay Area)

**PI (Physik Instrumente) L.P.**  
Sausalito, CA 94965  
www.pi-usa.us

### ITALY

**Physik Instrumente (PI) S. r. l.**  
Bresso  
www.pionline.it

### FRANCE

**PI France SAS**  
Aix-en-Provence  
www.pi.ws

### JAPAN

**PI Japan Co., Ltd.**  
Tokyo  
www.pi-japan.jp

### CHINA

**Physik Instrumente (PI Shanghai) Co., Ltd.**  
Shanghai  
www.pi-china.cn

### SOUTHEAST ASIA

**PI (Physik Instrumente) Singapore LLP**  
Singapore  
www.pi-singapore.sg  
For ID / MY / PH / SG / TH / VNM

### KOREA

**PI Korea Ltd.**  
Seoul  
www.pikorea.co.kr

### USA (West) & MEXICO

**PI (Physik Instrumente) L.P.**  
Irvine, CA 92620  
www.pi-usa.us

### UK & IRELAND

**PI (Physik Instrumente) Ltd.**  
Cranfield, Bedford  
www.physikinstrumente.co.uk

### NETHERLANDS

**PI Benelux B.V.**  
Sint-Oedenrode  
www.pi.ws/benelux

### SPAIN

**Micos Iberia S.L.**  
Vilanova i la Geltrú  
www.pimicos.es


### PI Japan Co., Ltd.

Osaka  
www.pi-japan.jp

**Physik Instrumente (PI Shanghai) Co., Ltd.**  
Beijing  
www.pi-china.cn

### TAIWAN

**Physik Instrumente (PI) Taiwan Ltd.**  
Taipei  
www.pi-taiwan.com.tw

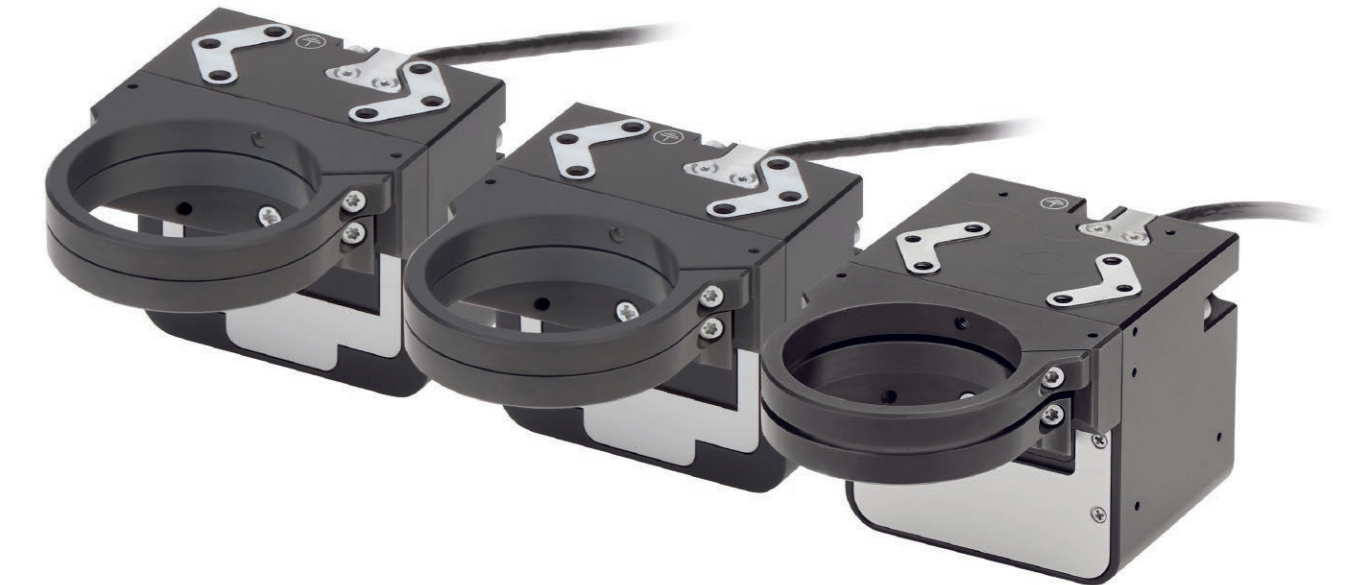

# P-725 PIFOC® PORTFOLIO

PIEZOBASED OBJECTIVE  
FOCUS DRIVE SERIES

BR077E P-725.xCDEZ PIFOC 06/2021 0.0. Subject to change without notice. © Physik Instrumente (PI) GmbH & Co. KG 2021

# Comprehensively Improved: PIFOC® Objective Focussing Drives

Large Travel Ranges, High Dynamics, Easy Handling

Be it in life sciences, materials microcopy, quality assurance, or semiconductor manufacturing – vertical positioning systems for objectives play an increasingly important role in numerous applications. Piezo-based, lever-amplified vertical axes have proven themselves as being the standard solution for this task for about 20 years now. PI's extensive portfolio of compact PIFOC® objective positioners has now undergone a complete redesign. This redesign wanted to achieve higher dynamics, travel ranges to 800 µm, and an improved usability. The result is an even more powerful range than before, which has the right solution for virtually all applications.

## Integration Made Easy


Simple Clamp Ring for Threaded Adapter

PI provides all the necessary accessories needed for the most common usage, which is screwing the PIFOC® into the revolving nosepiece. The aperture of the clamp ring is big enough for objectives up to M34. The threaded adapter for objectives of different sizes are fixed with Torx screws in the clamp ring. Furthermore, the PIFOC® can flexibly be integrated in platform assemblies. For this purpose, tolerated interfaces are placed in the base body.



# Main Benefits and Improved Usability

This is what Makes the Difference

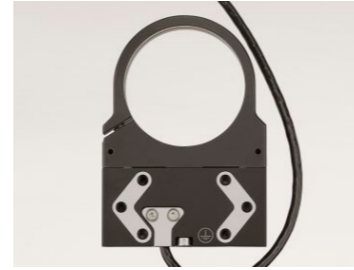

Highly flexible cables allow for tight bending radii / tolerated interfaces with threads for platform integration



New strain relief for robustness and longevity



Flat cross-sections of clamps for tight installation spaces



Matching spacer for same parfocal length / Slim cross-sections enable the use of directly adjacent lens mounts

### For Better Work Results

- Piezo drives with new lever designs
- Up to 800 µm travel range
- High productivity due to short data acquisition times/quick scanning
- High dynamics i.e. high resonant frequencies and short settling times
- Reliable performance based on precise motion and positioning

### Usability & Service

- Simple handling thanks to new strain relief and highly flexible cables
- Lowest possible restriction on the focus displacement in Z thanks to flat cross-sections of the clamps
- Always the right parfocal length when changing magnification thanks to adapted spacers
- No wear – nowhere: 24/7 operation thanks to wear-free, lever-amplified piezo drives
- Short delivery times
- Support by engineering and application experts



# An Unparalleled Wide Offer with Leading Key Features

From ≤ 18 µm to 7 mm Travel Range: The Right Solution for Any Application

	Travel range	Step&Settle (10% step, 150 g load)	Minimal incremental motion	Resonant frequency (@ 150 g)	Recommended electronics	Mounting interface(s)	Comments
 P-725.CDD	≤ 18 µm	5 ms	0,6 nm	450 Hz (@ 200 g)	E-709.1C1L E-754.1CD*	Objective thread adapter / Threads & location holes in base body	Shortest settling times
 P-726.1CD	100 µm	6 ms	1 nm	560 Hz (@ 200 g)	E-709.1C1L E-754.1CD*	Objective thread adapter	For highload objectives
 P-725.1CDE2	100 µm	14 ms	1 nm	290 Hz	E-709.1C1L E-754.1CD*	Objective thread adapter / Threads & location holes in base body	NEW: higher stiffness & improved usability
 P-725.4CDE2	400 µm	22 ms	4 nm	175 Hz	E-709.1C1L E-754.1CD*	Objective thread adapter / Threads & location holes in base body	NEW: higher stiffness & improved usability
 P-725.8CDE2	800 µm	39 ms	5 nm	110 Hz	E-709.1C1L E-754.1CD*	Objective thread adapter / Threads & location holes in base body	NEW: long travel
 ND722LAQ	2.000 µm	145 ms	5 nm	n.a.	E-861.1A1 (included in scope of delivery)	Objective thread adapter	Self-locking at rest due to NEXACT® technology
 V-308	7.000 µm	<15 ms (100nm; +-15nm)	10 nm	n.a.	C-414 G-910 ACS SPiiPlus + NPMpm (NanoPWM incl.)	Through-holes for mounting on flat surfaces	Adjustable weight force compensation to 1 kg

\* For special high performance requirements